

JC961 U.S. PTO

03-25-02

JC05 13134-PCT-US 22 MAR 2002

FORM PTO-1390 (REV. 9-2001)		U.S. DEPARTMENT OF COMMERCE PATENT AND TRADEMARK OFFICE	ATTORNEY'S DOCKET NUMBER <b>13134-PCT-US</b>
<b>TRANSMITTAL LETTER TO THE UNITED STATES DESIGNATED/ELECTED OFFICE (DO/EO/US) CONCERNING A FILING UNDER 35 U.S.C. 371</b>			U.S. APPLICATION NO. (if known, see 37 CFR 1.5) <b>10/089001</b>
INTERNATIONAL APPLICATION NO. PCT/GB00/03597	INTERNATIONAL FILING DATE 19 September 2000	PRIORITY DATE CLAIMED 24 September 1999	
TITLE OF INVENTION Assay			
APPLICANT(S) FOR DO/EO/US Orchid BioSciences, Inc.			
Applicant herewith submits to the United States Designated/Elected Office (DO/EO/US) the following items and other information:			
<p>1. <input checked="" type="checkbox"/> This is a <b>FIRST</b> submission of items concerning a filing under 35 U.S.C. 371.</p> <p>2. <input type="checkbox"/> This is a <b>SECOND</b> or <b>SUBSEQUENT</b> submission of items concerning a filing under 35 U.S.C. 371.</p> <p>3. <input type="checkbox"/> This is an express request to begin national examination procedures (35 U.S.C. 371(f)). The submission must include items (5), (6), (9) and (21) indicated below.</p> <p>4. <input type="checkbox"/> The US has been elected by the expiration of 19 months from the priority date (Article 31).</p> <p>5. <input checked="" type="checkbox"/> A copy of the International Application as filed (35 U.S.C. 371(c)(2))</p> <p>a. <input type="checkbox"/> is attached hereto (required only if not communicated by the International Bureau).</p> <p>b. <input checked="" type="checkbox"/> has been communicated by the International Bureau.</p> <p>c. <input type="checkbox"/> is not required, as the application was filed in the United States Receiving Office (RO/US).</p> <p>6. <input type="checkbox"/> An English language translation of the International Application as filed (35 U.S.C. 371(c)(2)).</p> <p>a. <input type="checkbox"/> is attached hereto.</p> <p>b. <input type="checkbox"/> has been previously submitted under 35 U.S.C. 154(d)(4).</p> <p>7. <input checked="" type="checkbox"/> Amendments to the claims of the International Application under PCT Article 19 (35 U.S.C. 371(c)(3))</p> <p>a. <input type="checkbox"/> are attached hereto (required only if not communicated by the International Bureau).</p> <p>b. <input type="checkbox"/> have been communicated by the International Bureau.</p> <p>c. <input checked="" type="checkbox"/> have not been made; however, the time limit for making such amendments has NOT expired.</p> <p>d. <input type="checkbox"/> have not been made and will not be made.</p> <p>8. <input type="checkbox"/> An English language translation of the amendments to the claims under PCT Article 19 (35 U.S.C. 371(c)(3)).</p> <p>9. <input type="checkbox"/> An oath or declaration of the inventor(s) (35 U.S.C. 371(c)(4)).</p> <p>10. <input type="checkbox"/> An English language translation of the annexes of the International Preliminary Examination Report under PCT Article 36 (35 U.S.C. 371(c)(5)).</p> <p><b>Items 11 to 20 below concern document(s) or information included:</b></p> <p>11. <input type="checkbox"/> An Information Disclosure Statement under 37 CFR 1.97 and 1.98.</p> <p>12. <input type="checkbox"/> An assignment document for recording. A separate cover sheet in compliance with 37 CFR 3.28 and 3.31 is included.</p> <p>13. <input type="checkbox"/> A <b>FIRST</b> preliminary amendment.</p> <p>14. <input type="checkbox"/> A <b>SECOND</b> or <b>SUBSEQUENT</b> preliminary amendment.</p> <p>15. <input type="checkbox"/> A substitute specification.</p> <p>16. <input type="checkbox"/> A change of power of attorney and/or address letter.</p> <p>17. <input type="checkbox"/> A computer-readable form of the sequence listing in accordance with PCT Rule 13ter.2 and 35 U.S.C. 1.821 - 1.825.</p> <p>18. <input type="checkbox"/> A second copy of the published international application under 35 U.S.C. 154(d)(4).</p> <p>19. <input type="checkbox"/> A second copy of the English language translation of the international application under 35 U.S.C. 154(d)(4).</p> <p>20. <input checked="" type="checkbox"/> Other items or information: <b>Partial International Search Report</b></p>			

JC13 Rec'd PCT/PTO 22 MAR 2002

U.S. APPLICATION NO. (if known, see 37 CFR 1.53) <b>10/089001</b>		INTERNATIONAL APPLICATION NO.		ATTORNEY'S DOCKET NUMBER	
--	--	-------------------------------	--	--------------------------	--

21. <input checked="" type="checkbox"/> The following fees are submitted: <b>BASIC NATIONAL FEE (37 CFR 1.492 (a) (1) - (5)):</b> Neither international preliminary examination fee (37 CFR 1.482) nor international search fee (37 CFR 1.445(a)(2)) paid to USPTO and International Search Report not prepared by the EPO or JPO ..... <b>\$1040.00</b>  International preliminary examination fee (37 CFR 1.482) not paid to USPTO but International Search Report prepared by the EPO or JPO ..... <b>\$890.00</b>  International preliminary examination fee (37 CFR 1.482) not paid to USPTO but international search fee (37 CFR 1.445(a)(2)) paid to USPTO ..... <b>\$740.00</b>  International preliminary examination fee (37 CFR 1.482) paid to USPTO but all claims did not satisfy provisions of PCT Article 33(1)-(4) ..... <b>\$710.00</b>  International preliminary examination fee (37 CFR 1.482) paid to USPTO and all claims satisfied provisions of PCT Article 33(1)-(4) ..... <b>\$100.00</b>  <b>ENTER APPROPRIATE BASIC FEE AMOUNT =</b>				<b>CALCULATIONS PTO USE ONLY</b>          <div style="border: 1px solid black; padding: 2px; width: fit-content; margin: 0 auto;">\$ 890.00</div>	
Surcharge of <b>\$130.00</b> for furnishing the oath or declaration later than <input type="checkbox"/> 20 <input type="checkbox"/> 30 months from the earliest claimed priority date (37 CFR 1.492(e)).				\$	
CLAIMS	NUMBER FILED	NUMBER EXTRA	RATE	\$	
Total claims	30 - 20 =	10	x \$18.00	\$ 180	
Independent claims	10 - 3 =	7	x \$84.00	\$ 588	
MULTIPLE DEPENDENT CLAIM(S) (if applicable)				+ \$280.00	
<b>TOTAL OF ABOVE CALCULATIONS =</b>				\$	
<input type="checkbox"/> Applicant claims small entity status. See 37 CFR 1.27. The fees indicated above are reduced by 1/2.				\$	
<b>SUBTOTAL =</b>				\$	
Processing fee of <b>\$130.00</b> for furnishing the English translation later than <input type="checkbox"/> 20 <input type="checkbox"/> 30 months from the earliest claimed priority date (37 CFR 1.492(f)).				\$	
<b>TOTAL NATIONAL FEE =</b>				\$	
Fee for recording the enclosed assignment (37 CFR 1.21(h)). The assignment must be accompanied by an appropriate cover sheet (37 CFR 3.28, 3.31). <b>\$40.00</b> per property +				\$	
<b>TOTAL FEES ENCLOSED =</b>				\$ 1658	
				Amount to be refunded:	\$
				charged:	\$

a. ☐ A check in the amount of \$ \_\_\_\_\_ to cover the above fees is enclosed.

b. ☒ Please charge my Deposit Account No. 11-0171 in the amount of \$ 1658.00 to cover the above fees.  
A duplicate copy of this sheet is enclosed.

c. ☒ The Commissioner is hereby authorized to charge any additional fees which may be required, or credit any  
overpayment to Deposit Account No. 11-0171. A duplicate copy of this sheet is enclosed.

d. ☐ Fees are to be charged to a credit card. **WARNING:** Information on this form may become public. **Credit card  
information should not be included on this form.** Provide credit card information and authorization on PTO-2038.

**NOTE:** Where an appropriate time limit under 37 CFR 1.494 or 1.495 has not been met, a petition to revive (37 CFR  
1.137 (a) or (b)) must be filed and granted to restore the application to pending status.

SEND ALL CORRESPONDENCE TO:

Kalow & Springut LLP  
488 MADISON AVE  
New York NY 10022

*William D. Schmidt*  
SIGNATURE  
William D. Schmidt  
NAME  
39,492  
REGISTRATION NUMBER

Rec'd PCT/PTO 06 JAN 2003

10089001.022003

#4

25840.1  
PATENT

13134 PCT US

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Applicant: Bayliffe, et al.

Serial No.: 10/089,001

International Application No.: PCT/GB00/03597

Date Filed: March 22, 2002

International Filing Date: 19 September 2000

International Priority Date: 24 September, 2000

Title: "ASSAY"

Customer No.: 23719

PATENT TRADEMARK OFFICE

Kalow & Springut LLP  
488 Madison Avenue, 19th Floor  
New York, New York 10022  
(212) 813-1600

January 6, 2003

Commissioner for Patents  
Washington, D.C. 20231

TRANSMITTAL OF SEQUENCE LISTING

Sir:

This is a reply to the Notice to Comply With Requirements for Patent Applications Containing Nucleotide Sequence and/or Amino Acid Sequence Disclosures issued in connection with the above-identified patent application on June 4, 2002. The deadline for reply to this Notice is January 6, 2003. Accordingly, this reply is timely filed. Enclosed also please find a petition for extension of time.

Certificate of Mailing Under 37 C.F.R. 1.8

I hereby declare that this correspondence is being deposited with the United States Postal Service as Express Mail Label No.: **EV 035748108 US** in an envelope addressed to: Commissioner for Patents, Box PCT, Washington, D.C. 20231.

*Jan 6-2003*  
January 6, 2003

*Sylvia Gonzalez*  
Sylvia Gonzalez

Applicants: Bayliffe et al  
International Application No.: PCT/BG00/03597  
International Filing Date: 19 September 2000  
Page -2- (Notification of Missing Requirements under 35 U.S.C. 371)  
Sequence Listing

13134-PCT-US

Please amend the above-identified application as follows:

**IN THE SPECIFICATION:**

Please remove pages 1-22 of the sequence listing and replace them with pages 1-26 of the attached substitute sequence listing.


**REMARKS**

This amendment replaces the original sequence listing with the attached substitute sequence listing. A computer readable copy of the sequence is enclosed. A paper copy is also being submitted. The computer readable copy of the substitute sequence listing and the attached paper copy are identical and meet the requirements of 37 CFR 1.821-1.825. No new matter has been added.

If any additional fees are due, or an overpayment has been made, please charge, or credit, our Deposit Account No. 11-0171 for such sum.

If the Examiner has any questions regarding the present application, the Examiner is cordially invited to contact Applicants' attorney at the telephone number provided below.

Respectfully submitted,

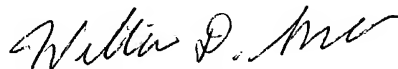

William D. Schmidt  
Registration No.: 39,492  
Attorney for Applicant

Kalow & Springut LLP  
(212) 813-1600



4324.1

13134-PCT-US

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Applicants: Bayliffe et al. Examiner: To be assigned  
Serial No.: To be assigned Art Unit: To be assigned  
Filed: Herewith  
Title: Assay

Kalow & Springut LLP  
488 Madison Avenue  
New York, New York 10022

March 22, 2002

Commissioner for Patents  
Washington, D.C. 20231

PRELIMINARY AMENDMENT

Dear Sir:

Prior to examination on the merits of this national stage application under the Patent Cooperation Treaty, please enter the following amendments in the above-identified application.

AMENDMENT

IN THE SPECIFICATION:

On page 1 after the title, please insert the following paragraph:

Certificate of Mailing Under 37 C.F.R. 1.10

I hereby declare that this correspondence is being deposited with the United States Postal Service via Express Mail Label No. E.V03577643545 in an envelope addressed to: Assistant Commissioner for Patents, Washington, D.C.

Date: 3/22/02 Name: J. Colwell

Applicant: Bayliffe et al.  
Serial No.: Not yet assigned  
Docket 13134-PCT-US  
Preliminary Amendment – March 22, 2002  
Page 2 of 3

**--Cross-Reference to Related Applications**

This application is a national stage filing under the Patent Cooperation Treaty (PCT) for PCT international application number PCT/GB00/03597 filed on 19 September 2000, published under PCT Article 21(2) in English as WO 01/21833 on 29 March 2001. Both the international application and this national stage application claim the benefit of the filing date of GB 9922527.8 filed 24 September 1999.--

**IN THE CLAIMS**

Please amend claims 6-9 to read as follows:

6. (Amended) A method of claim 1 which further comprises the use of one or more common amplification primer(s) in the presence of appropriate nucleotide triphosphates and an agent for polymerization, and subjecting the mixture to PCR amplification such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a corresponding PCR amplification product.
7. (Amended) A method of claim 2 wherein two or more diagnostic primers are used as a multiplex.
8. (Amended) A method of claim 2 wherein all of the diagnostic primers are used in a single multiplex reaction.
9. (Amended) A method of claim 1 which further comprises the use of one or more control primers.

Applicant: Bayliffe et al.  
Serial No.: Not yet assigned  
Docket 13134-PCT-US  
Preliminary Amendment – March 22, 2002  
Page 3 of 3

**REMARKS**


Entry of this preliminary amendment is respectfully requested. Pursuant to 37 CFR §1.121, a marked up version of the amendment to the specification and to the amended claims is submitted herewith.

The claim amendments reflect changes to conform the claims to U.S. practice. The amendments do not add new matter.

This is a PCT national stage filing under 35 USC §371, claiming a right of priority to an earlier filed foreign application GB 9922527.8 filed 24 September 1999, in accordance with 35 USC §363 and 35 USC §365. The amendment to the specification identifies the present application as a national stage filing in the U.S. under the PCT (for PCT/GB00/03597 filed 19 September 2000) and also claims priority to GB 9922527.8 filed 24 September 1999.

Applicants submit no fee is required in connection with the filing of this Preliminary Amendment. If any fee is deemed necessary, authorization is hereby given to charge the amount of any such fee to Deposit Account No. 11-0171.

Respectfully submitted,

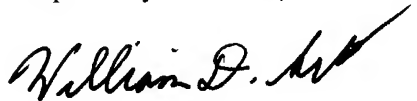

William D. Schmidt  
Attorney for Applicant(s)  
Registration No.: 39,492

4379.1

13134-PCT-US

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

Applicants: Bayliffe et al. Examiner: To be assigned  
Serial No.: To be assigned Art Unit: To be assigned  
Filed: Herewith  
Title: Assay

Kalow & Springut LLP  
488 Madison Avenue  
New York, New York 10022

March 22, 2002

Commissioner for Patents  
Washington, D.C. 20231

MARKED UP CLAIMS

Dear Sir:

In accordance with 37 CFR §1.121 (c) the following marked up claims are submitted to accompany the amendment filed concurrently for the application identified above.

6. (Amended) A method [as claim in any one of the previous claims and] of claim 1 which further comprises the use of one or more common amplification primer(s) in the presence of appropriate nucleotide triphosphates and an agent for polymerization, and subjecting the mixture to PCR amplification such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a corresponding PCR amplification product.

Certificate of Mailing Under 37 C.F.R. 1.10

I hereby declare that this correspondence is being deposited with the United States Postal Service via Express Mail Label No. E103576432US in an envelope addressed to: Assistant Commissioner for Patents, Washington, D.C.

Date: 3/22/02 Name: J. Kalow

7. (Amended) A method [as claimed in any one of claims] of claim 2[-6 and] wherein two or more diagnostic primers are used as a multiplex.
8. (Amended) A method [as claimed in any one of claims] of claim 2[-6 and] wherein all of the diagnostic primers are used in a single multiplex reaction.
9. (Amended) A method [as claimed in any one of the previous claims and] of claim 1 which further comprises the use of one or more control primers.

Respectfully submitted,



William D. Schmidt  
Attorney for Applicant(s)  
Registration No. 39,492

Kalow & Springut LLP  
(212) 813-1600

SEQUENCE LISTING

<110> Bayliffe, Andrew  
 Docter, Eelco  
 Kelly, Stephen James  
 Robertson, Nancy Hastings

<120> Assay

<130> 13134PCTUS

<140> US 10/089,001

<141> 2002-03-22

<150> PCT/GB00/03597

<151> 2002-09-19

<150> GB 9922527.8

<151> 1999-09-24

<160> 108

<170> FastSEQ for Windows Version 4.0

<210> 1

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature

<222> (1)...(26)

<223> n= A, C, G, or T, wherein 0-20 n may be absent

<400> 1

nnnnnnnnnn nnnnnnnnnn nnnnnntggt aaagac

36

<210> 2

<211> 66

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature

<222> (1)...(56)

<223> n = A, C, G, or T, wherein 0-40 n may be absent

<400> 2

nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnttaa 60

aaaagc

66

<210> 3  
 <211> 66  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(56)  
 <223> n = A, C, G, or T, wherein 0-40 n may be absent

<400> 3  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnaaaa 60  
 aaaatc 66

<210> 4  
 <211> 26  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 4  
 taattcccca aatccctgtt aaagac 26

<210> 5  
 <211> 88  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(60)  
 <223> n = A, C, G, or T, wherein 0-55 n may be absent

<400> 5  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn 60  
 taattcccca aatccctgtt aaaaaagc 88

<210> 6  
 <211> 90  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(60)  
 <223> n = A, C, G, or T, wherein 0-55 n may be absent

<400> 6  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn 60  
 taattcccca aatccctggt aaaaaaatc 90

<210> 7  
 <211> 53  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 7  
 gttaatcatt cagctactac gcacctaatt ccccaaattc ctgttaaaaa agc 53

<210> 8  
 <211> 89  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 8  
 gactgtacga tactcattta tatgaagtca gctacttacc tatagaacgc ttgctagttt 60  
 aattcccca atccctgtta aaaaaaatc 89

<210> 9  
 <211> 35  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 9  
 gtacataaaa caagcatcta ttgaaaatat ctgac 35

<210> 10  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 10  
 nnnnnnnnnn nnnnnnnnnn nnnnnngcaa cagtca 36

<210> 11  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence



<220>  
<223> Primer

<221> misc\_feature  
<222> (1)...(26)  
<223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 11  
nnnnnnnnnn nnnnnnnnnn nnnnnnttgg taatta 36

<210> 12  
<211> 36  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<221> misc\_feature  
<222> (1)...(26)  
<223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 12  
nnnnnnnnnn nnnnnnnnnn nnnnnnatag ttctct 36

<210> 13  
<211> 36  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<221> misc\_feature  
<222> (1)...(26)  
<223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 13  
nnnnnnnnnn nnnnnnnnnn nnnnnnggga tccatc 36

<210> 14  
<211> 36  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<221> misc\_feature  
<222> (1)...(26)  
<223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 14  
nnnnnnnnnn nnnnnnnnnn nnnnnnaaac accatt 36

<210> 15

<211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 15  
 nnnnnnnnnn nnnnnnnnnn nnnnnngtct tacgca 36

<210> 16  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 16  
 tcttgggatt caataacttt gcaacagtca 30

<210> 17  
 <211> 27  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 17  
 tctcgaattt tctatTTTTg gtaatta 27

<210> 18  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 18  
 agtttgcaga gaaagacaat atagttctct 30

<210> 19  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 19  
 tgatcactcc actgttcata gggatccatc 30

<210> 20  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 20  
 gtatctatat tcatcatagg aaacaccatt 30

<210> 21  
 <211> 29  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 21  
 gaacatttcc tttcagggtg tcttacgca 29

<210> 22  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 22  
 nnnnnnnnnn nnnnnnnnnn nnnnnnaaac accaca 36

<210> 23  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(6)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 23  
 nnnnnnnnnn nnnnnnnnnn nnnnnngcga tagact 36

<210> 24  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 24  
 nnnnnnnnnn nnnnnnnnnn nnnnnngaag tattga 36  
  
 <210> 25  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 25  
 nnnnnnnnnn nnnnnnnnnn nnnnnnatca tcctgt 36  
  
 <210> 26  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 26  
 nnnnnnnnnn nnnnnnnnnn nnnnnntctg tgagtt 36  
  
 <210> 27  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 27  
 nnnnnnnnnn nnnnnnnnnn nnnnnnttct tgctga 36  
  
 <210> 28

<211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 28  
 nnnnnnnnnn nnnnnnnnnn nnnnnngctc gttgtt 36

<210> 29  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 29  
 agcctatgcc tagataaatc gcgatagact 30

<210> 30  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 30  
 tgccatgggg cctgtgcaag gaagtattga 30

<210> 31  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 31  
 cctatgcact aatcaaagga atcatcctgt 30

<210> 32  
 <211> 30  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 32  
 tatttttatt tcagatgcga tctgtgagtt 30

<210> 33  
 <211> 31  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 33  
 ttattcacct tgctaaagaa attcttgctg a 31

<210> 34  
 <211> 25  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 34  
 gctaaagaaa ttcttgctcg ttgtt 25

<210> 35  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 35  
 nnnnnnnnnn nnnnnnnnnn nnnnnnagtt gttgta 36

<210> 36  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 36  
 nnnnnnnnnn nnnnnnnnnn nnnnnnaggg ttctctg 36

<210> 37  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 37  
 nnnnnnnnnn nnnnnnnnnn nnnnnnttgt tctacc 36  
  
 <210> 38  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 38  
 nnnnnnnnnn nnnnnnnnnn nnnnnngaaa ataact 36  
  
 <210> 39  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 39  
 nnnnnnnnnn nnnnnnnnnn nnnnnntaaa cctaga 36  
  
 <210> 40  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 40  
 nnnnnnnnnn nnnnnnnnnn nnnnnnaaaa gatagc 36  
  
 <210> 41

<211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 41  
 nnnnnnnnnn nnnnnnnnnn nnnnnncagg gaagca

36

<210> 42  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 42  
 nnnnnnnnnn nnnnnnnnnn nnnnnnaagc cagtta

36

<210> 43  
 <211> 29  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 43  
 ttcaagatag aaagaggaca gttgttgta

29

<210> 44  
 <211> 23  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 44  
 ccttcttctt ctcagggttc ctg

23

<210> 45  
 <211> 27  
 <212> DNA  
 <213> Artificial Sequence

<220>



<223> Primer

<400> 45

caccatctca ttctgcattg ttctacc

27

<210> 46

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 46

gcctggcacc attaaagaaa ataact

26

<210> 47

<211> 24

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 47

atgccaacag aaggtaaacc taga

24

<210> 48

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 48

caaactctcc agtctgttta aaagatagc

29

<210> 49

<211> 23

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 49

ggaagaactg gatcagggaa gca

23

<210> 50

<211> 26

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 50  
 ttaggatttt tctttgaagc cagtta 26  
  
 <210> 51  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 51  
 nnnnnnnnnn nnnnnnnnnn nnnnnntggt ctacga 36  
  
 <210> 52  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 52  
 nnnnnnnnnn nnnnnnnnnn nnnnnntatt tagtga 36  
  
 <210> 53  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (2)...(26)  
 <223> n = A, C, G, or T, wherein 0-20 n may be absent  
  
 <400> 53  
 nnnnnnnnnn nnnnnnnnnn nnnnnncaca ctgacg 36  
  
 <210> 54  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature

<222> (1)...(26)

<223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 54

nnnnnnnnnn nnnnnnnnnn nnnnnnctgc caaata

36

<210> 55

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature

<222> (1)...(26)

<223> n = A, C, G, or T, wherein 0-20 n may be absent

<400> 55

nnnnnnnnnn nnnnnnnnnn nnnnnntatc cacttg

36

<210> 56

<211> 50

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 56

tagccattga tgacggagcg atgttttttc tggagattta tgttctacga

50

<210> 57

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 57

gatttatgtt ctatggaatc tttttatatt tagtga

36

<210> 58

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 58

tgagagaaggt ggaatcacac tgacg

25

<210> 59

<211> 28

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 59

aagctctgaa ttacatact gccaaata

28

<210> 60

<211> 33

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 60

aaagatgata agacttacca agctatccac ttg

33

<210> 61

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 61

taatctctac caaatctgga tactatacc

29

<210> 62

<211> 29

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 62

taatctctac caaatctgga tactatacc

29

<210> 63

<211> 25

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 63

gaattcccaa acttttagag acatc

25

<210> 64

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 64

cttgatggta agtacatggg tttttcttat

30

<210> 65

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 65

ccagacttca cttctaatta tgattatggg

30

<210> 66

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 66

ttgtggatca aatttcagtt gacttgatc c

31

<210> 67

<211> 30

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 67

gacttcactt ctaatgatga ttatgggaga

30

<210> 68

<211> 37

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 68

gtttcacata gtgtatgacc ctctatatac actcatt

37

<210> 69

<211> 31

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<400> 69  
atctaaaatt ggagcaatgt tgtttttgac c 31

<210> 70  
<211> 31  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<400> 70  
atctaaaatt ggagcaatgt tgtttttgac c 31

<210> 71  
<211> 37  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<400> 71  
gtttcacata gtgtatgacc ctctatatac actcatt 37

<210> 72  
<211> 29  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<400> 72  
ttttgctgtg agatctttga cagtcattt 29

<210> 73  
<211> 30  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<400> 73  
tttgtttatt gctccaagag agtcatacca 30

<210> 74  
<211> 49  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<400> 74

gactgactga ctgactgaaa tggagacttt ttgtttatgt ggttactaa

49

<210> 75  
 <211> 35  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 75  
 gtatgataga gattatatgc aataaaacat taaca

35

<210> 76  
 <211> 34  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 76  
 tgtgtctaatt attgattcta ctgtacaata ataa

34

<210> 77  
 <211> 35  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 77  
 atttttccaa acttcattag aactgatcta ttgac

35

<210> 78  
 <211> 27  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 78  
 cacagtagct tacccataga ggaaaca

27

<210> 79  
 <211> 35  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 79  
 atttttccaa acttcattag aactgatcta ttgac

35

<210> 80  
 <211> 23  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 80  
 gctcacctgt ggtatcactc caa 23

<210> 81  
 <211> 34  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 81  
 aatcaaacta tgttaaggga aataggacaa ctaa 34

<210> 82  
 <211> 39  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 82  
 cgattcgatt cagttttctg tggtttctta gtgtttgga 39

<210> 83  
 <211> 33  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 83  
 gtaatttttt tacatgaatg acatttacag caa 33

<210> 84  
 <211> 33  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 84  
 ggaaattatt tgtttaacaa taaaacaatg gaa 33

<210> 85  
 <211> 26



<212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<400> 85  
 ccaacaacac ctccaatacc agtaac

26

<210> 86  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 86  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn aaagac

36

<210> 87  
 <211> 96  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(90)  
 <223> n = A, C, G, or T wherein 0-75 n may be absent

<400> 87  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn 60  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn aaaagc 96

<210> 88  
 <211> 96  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(90)  
 <223> n = A, C, G, or T wherein 0-75 n may be absent

<400> 88  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn 60  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn aaaatc 96

<210> 89

<211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 89  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn accatt

36

<210> 90  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 90  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn tacgca

36

<210> 91  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 91  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn tccatc

36

<210> 92  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer

<221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 92

nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn taatta

36

<210> 93  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
 <220>  
 <223> Primer  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 93  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn cagtca

36

<210> 94  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
 <220>  
 <223> Primer  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 94  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn ttctct

36

<210> 95  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
 <220>  
 <223> Primer  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 95  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn accaca

36

<210> 96  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
 <220>  
 <223> Primer  
 <221> misc\_feature  
 <222> (1)...(30)

<223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 96

nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn gttgta

36

<210> 97

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature

<222> (1)...(30)

<223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 97

nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn gatagc

36

<210> 98

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature

<222> (1)...(30)

<223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 98

nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn cctaga

36

<210> 99

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature

<222> (1)...(30)

<223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 99

nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn ataact

36

<210> 100

<211> 36

<212> DNA

<213> Artificial Sequence

<220>

<223> Primer

<221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 100  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn ttcctg 36  
  
 <210> 101  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 101  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn tctacc 36  
  
 <210> 102  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 102  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn gaagca 36  
  
 <210> 103  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 103  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn cagtta 36  
  
 <210> 104  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence

<220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 104  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn ctacga 36  
  
 <210> 105  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 105  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn tagtga 36  
  
 <210> 106  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 106  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn ctgacg 36  
  
 <210> 107  
 <211> 36  
 <212> DNA  
 <213> Artificial Sequence  
  
 <220>  
 <223> Primer  
  
 <221> misc\_feature  
 <222> (1)...(30)  
 <223> n = A, C, G, or T wherein 0-20 n may be absent  
  
 <400> 107  
 nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn caaata 36  
  
 <210> 108

<211> 36  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Primer

<221> misc\_feature  
<222> (1)...(30)  
<223> n = A, C, G, or T wherein 0-20 n may be absent

<400> 108  
nnnnnnnnnn nnnnnnnnnn nnnnnnnnnn cacttg

36

WO 01/21833

PCT/GB00/03597

- 1 -

ASSAY

5 This invention relates to a diagnostic method for the simultaneous detection of the variable length polythymidine (polyT) tract alleles in the CFTR gene that are associated with the phenotypic modulation of selected CFTR mutations (Friedman KJ. et al. Hum Mutat. 10: 108-115 (1997)).

10 The invention also relates to a diagnostic method further comprising the simultaneous detection of one or more of the DF508 (Kerem B. et al. Science 245: 1073-1080 (1989)), 621+1G>T (Rozen R. et al. Am. J. Med. Genet. 42: 360-364 (1992)), G542X, 3659delC, A455E, DI507 (Kerem B. et al. Proc. Nat. Acad. Sci. 87: 8447-8451 (1990)), 3849+10kbC>T (Highsmith WE. et al. New Eng. J. Med. 331: 974-980 (1994)), N1303K (Osborne L. et al. Am. J. Hum. Genet. 48: 608-612 (1991)), 1717-1G>A (Kerem B. et al. Proc. Nat. Acad. Sci. 87: 8447-8451 (1990)), Guillermit H. et al. Hum. Genet. 85: 450-453 (1990)), 1078delT (Claustres M. et al. Genomics 13: 907-908 (1992)), W1282X (Vidaud M et al. Hum. Genet. 15 85: 446-449 (1990)), R347P, R117H (Dean M. et al. Cell 61: 863-870 (1990)), G551D, R553X (Cutting GR. et al. New Eng. J. Med. 323: 1685-1689 (1990)), S1251N (Gasparini P. et al. Hum. Mutat. 2: 389-394 (1993)), R1162X, R334W (Gasparini P. et al. Genomics 10: 193-200 (1991)), 2183AA>G, (Bozon D. et al. Hum. Mutat. 3: 330-332 (1994)), and the E60X (Will K. et al. Hum. Mutat. 5: 210-20 (1995)) mutations in the human cystic fibrosis conductance regulator (CFTR) gene that are prevalent in populations of European descent using the amplification refractory mutation system, (ARMS).

25 The invention also relates to a diagnostic method further comprising the simultaneous detection of one or more of the G85E (Chalkley G. et al. J. Med. Genet. 28: 875-877 (1991)), 405+1G>A (Dork T. et al. Hum Mol Genet. 2: 1965-1966 1993)), S549R (Kerem B. et al. Proc. Nat. Acad. Sci. 87: 8447-8451 (1990)), (or, depending on which of the two mutations produce this phenotype is tested for, Sangiuolo. et al. Genomics 9: 788-789 (1991)), W1089X (Shoshani T. et al. Hum. Molec. Genet. 3: 657-658 (1994)), D1152H (Feldmann D. et al. Clin. Chem. 41: 1675 (1995)) mutations in the human CFTR gene that are prevalent in populations of Middle Eastern descent using ARMS.

30 The invention also relates to mutation specific and allele specific primers for use in the method, to diagnostic kits containing these primers and to techniques for detecting primer specific amplification products.



Cystic fibrosis transmembrane conductance regulator (CFTR) functions as a chloride channel and controls the regulation of other transport pathways. The CFTR gene maps to chromosome 7q (Riordan J R. et al. Science 245: 1066-1073 (1989)). Mutations in the CFTR gene cause cystic fibrosis (CF) formerly known as cystic fibrosis of the pancreas. Cystic  
5 fibrosis is an autosomal recessive disorder mainly among Caucasians with a frequency of approximately 1/2 500 (Welsh MJ. et al. in Scriver et al. (eds): The Metabolic Basis of Inherited Disease, Vol 3: 7th ed. McGraw-Hill; New York, pp3799-3876 (1995)). Cystic fibrosis disrupts exocrine function of the pancreas, intestinal glands (meconium ileus), biliary tree (biliary cirrhosis), bronchial glands (chronic bronchopulmonary infection with  
10 emphysema), and sweat glands (high sweat electrolyte with depletion in a hot environment). Infertility occurs in both males and females. Many of the clinical and pathological findings are thought to be attributable to a generalised defect in mucus secretion due to an abnormality in the chloride channel.

The treatment of CF is complex, costly and time-consuming. Most children with CF in  
15 the developed world are followed in CF clinics which are typically multidisciplinary, involving physicians, nurses, social workers, nutritionists and physiotherapists. Treatment is directed at improving nutrition through the use of replacement pancreatic enzymes and vitamins, as well as a high-energy, high-protein and liberal-fat diet. For those not responding to this approach, enteral supplementation by nightly nasogastric, gastrostomy or jejunostomy  
20 infusion of high-energy diets has been used. Pulmonary treatment includes antibiotic therapy, either maintained continuously or reserved for exacerbations, and chest physiotherapy consisting of postural drainage, percussion, vibration and assisted coughing.

Approximately 70% of the mutations in CF patients correspond to a specific deletion of 3 basepairs, which results in the loss of a phenylalanine residue at amino acid position 508  
25 of the CFTR protein (Kerem B. et al. Am. J. Hum. Genet. 44: 827-834 (1989)).

Mutations of the CFTR gene are believed to produce an abnormal protein as a component of the chloride channel gate at the cell surface and are classified according to their phenotypic manifestations. Class I are nonsense mutations, resulting in the introduction of a stop codon, and in no synthesis of CFTR e.g. G542X, R553X, W1282X. Class II mutations  
30 cause a block in CFTR assembly in the endoplasmic reticulum or affect transport to the cell membrane, these include  $\Delta F508$ , A455E, and P574H. Class III and IV mutations give rise to defective chloride channel activity or regulation e.g. R117H, which alters amino acid residues

in the first CFTR transmembrane domain. Class V mutations result in modulations in CFTR synthesis. In cases of CF involving class I and II CFTR mutations (most commonly  $\Delta F508$  homozygotes) chronic pulmonary disease is observed at variable levels with chronic pancreatitis. However, some CF patients present with "mild" symptoms ranging from

5 sinopulmonary problems with pancreas sufficiency (PS-CF), to congenital bilateral absence of vas deferens (CBAVD). Mutations associated with PS-CF include R334W and R117H; the latter mutation accounts for approximately 0.8% of mutant alleles in Caucasian CF patients (Tsui LC. Trends Genet. 57: 392-398 (1992)).

There is a highly conserved polyT at the end of intron 8 of the CFTR gene. This polyT

10 tract plays an important role in aiding splice branch site recognition, and identifying and using the splice acceptor site (Smith CW. et al. Nature. 342: 243-247 (1989)). There are three allelic variants depending on the number of thymidines (5, 7, or 9) with allele frequencies in the general population of 5%, 85%, 10%, and respectively (Chu C-S. et al. EMBO J. 10: 1355-1363 (1991)). The number of thymidines determines the efficiency by which the intron 8

15 splice acceptor site is used. The efficiency decreases when a shorter stretch of thymidine residues is found. The T5 allele results in the most inefficient use of this splice acceptor site (Teng, H. et al. Hum. Mol. Genet. 6: 85-90 (1997)) and so CFTR transcripts from a T5 allele will lack exon 9 sequence. Individuals homozygous for the 5T allele were found to generate 90% CFTR mRNA with exon9 skipped (exon9<sup>-</sup>). While 7T homozygotes generated <25%

20 exon9<sup>-</sup> transcripts, and 9T <15% (Lissens W. et al. Hum. Reprod. 11: Supp4 55-77 (1996)).

There is an association between phenotype and the particular polyT allelic background for some CF mutations, for example, if a CFTR gene with the R117H mutation harbors a T5 allele, the mutant gene will be responsible for CF. An R117H mutant CFTR gene that harbors a T7 allele can either result in CF or CBAVD (Kiesewetter S. et al. Nature Genet. 5: 274-278,

25 (1993)).

The existence of polyT splice site variants may therefore explain the clinical heterogeneity observed in some mild cases of CF. It is also possible that the CF clinical phenotypes associated with mutations other than  $\Delta F508$  result from both the mutation in the CFTR gene and the sequence at exon9 splice site acceptor site.

30 Most CFTR mutations occur on a particular polyT background, most commonly the 9T allele;  $\Delta F508$  is in linkage disequilibrium with the 9T allele and the mild CF mutation R117H is always found on a 5T or 7T background. The mutation R117H has been shown to

owe its phenotypic heterogeneity (PS CF or CBAVD) to the presence of the 5T or 7T polyT alleles (Friedman KJ. et al. Hum Mutat. 10: 108-115 (1997)). Compound heterozygotes having the genotype  $\Delta F508/R117H$  may present as PS-CF, CBAVD or be asymptomatic. Since the mutation  $\Delta F508$  is in linkage disequilibrium with the 9T allele, it has been

5 demonstrated that the R117H phenotype contrasts with respect to its associated polyT haplotype (R117H-5T or R117-7T). The incidence of 5T alleles is also associated with other clinical phenotypes, for example, disseminated bronchiectasis (Pignatti PF. et al. Am. J. Hum. Genet. 58: 889-892 (1996)). Hence the human CFTR gene intron 8 5T, 7T and 9T allele

10 discriminatory test described and disclosed here is a test that is particularly useful when used in combination with any other cystic fibrosis diagnostic test.

Several cohort studies of screened and unscreened subjects have suggested that people identified with CF in the presymptomatic phase do better than those in whom a diagnosis is made because of symptoms. For example, one study (Dankert-Roelse JE. et al. J. Pediatr. 114: 362-367 (1989)) showed that, 88% of screened children but only 60% of unscreened children

15 were still alive at age 11 years.

Diagnostic techniques for phenotypic manifestations of CF include measuring sweat chloride. Elevated sweat chloride concentrations (greater than 60 mmol/l) are almost exclusively observed in patients with CF. However, it is difficult to obtain sufficient sweat from newborns, even after stimulating localised sweating by administering pilocarpine into

20 the skin it is often impossible to collect sufficient sweat for accurate analysis. Sweat chloride concentrations can be measured, for example, by using the Lazar ISM-146 Micro Chloride electrode. Measurement of elevated immunoreactive trypsinogen (IRT) in a dried blood spot is another screening method for CF. False positives and false negatives are known to occur, with false negatives occurring more frequently in neonates with meconium ileus. The positive

25 predictive value of the test is only 1-7% (Ryley HC. J. Clin. Pathol. 41: 726-729 (1988), Edminson PD. et al. Scand. J. Gastroenterol. Suppl. 143: 13-18 (1988)). A false positive rate of 93-99% could and likely does generate considerable anxiety, which may be long-lasting (Ryley HC. J. Clin. Pathol. 41: 726-729 (1988)). Furthermore, elevations of IRT decline after the first few months of life, so while exact timing of specimen collection in the neonatal

30 period is not critical, the collection of a second screening specimen to follow-up an initial abnormal screen should occur no earlier than 21 days, to avoid an increased number of false positives, and no later than 60 days, to reduce the risk of false negatives. Some newborn

screening programs rely on other tests. The "BM meconium test" identifies infants with high albumin content in the stool resulting from pancreatic insufficiency but it has very low sensitivity (Naylor EW. Semin. Perinatol. 9: 232-249 (1985)).

5 Friedman et al. describe an allele-specific amplification polyT assay (Friedman KJ. et al. Hum. Mutat. 10: 108-115 (1997)) designed to screen individuals for CFTR intron 8 5T, 7T, 9T alleles and heterozygotic polyT allelic combinations. However, the assay described by Friedman requires the use of each allele-specific primer in a separate reaction.

10 In our European Patent No. 0 332 435 B1 we disclose the Amplification Refractory Mutation System (ARMS). This simple and elegant method permits the detection of point mutations via allele-specific amplification of target sequences. In EP-0 332 435 we disclose and claim the application of ARMS to a variety of inherited and/or acquired genetic disorders.

Multiplex allele-specific amplification assays have been described for the simultaneous analysis of several CF causing mutations (for example, see Ferrie RM. et al. 15 Am. J. Hum. Genet. 51: 251-262 (1992) and Robertson NH. et al. Eur. Respir. J. 12: 477-482 (1998)) as disclosed in our European Patent No. 0 928 832. The combination of different primers for the simultaneous detection of two or more point mutations is termed "multiplexing". However, the design of robust and accurate allele-specific multiplex tests is not straightforward. A person skilled in the art would still need to contend with inter- and 20 intra- primer reactions when several alleles are amplified simultaneously in one reaction; there may be further interactions between amplification products for one allele with the primers for another allele. Furthermore, there will be inequalities between the efficiencies amplification products from different alleles.

25 We now provide allele-specific primers for the 5T, 7T and 9T alleles in intron 8 of the CFTR gene that are designed for use in combination in a multiplex assay and function by bringing about an allele-specific size differential between CFTR intron 8 5T, 7T and 9T derived amplification products.

Therefore according to a first aspect of the invention we now provide a diagnostic method for the detection of the 5T, 7T and 9T alleles in intron 8 of the human CFTR gene 30 which method comprises contacting a test sample of nucleic acid from an individual with a multiplex of diagnostic primers comprising (i) 5T variant primer 5'(N)nAAAGAC3', (ii) 7T variant primer 5'(N\*)n\*(N)nAAAAGC3' and (iii) 9T variant primer 5'(N\*)n\*(N)nAAAATC3', wherein N represents additional nucleotides which base pair with

the corresponding genomic sequence in the respective allele and n is an integer between 10 and 30 and N\* represents additional non-homologous nucleotides which do not base pair with the corresponding genomic sequence in the respective allele and n\* is an integer between 5 and 60, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

It will be appreciated that N\* may be any other molecule that serves to reduce the electrophoretic mobility of an extension product such as 3'hexaethylene glycol (HEG) or a combination of such molecules and additional non-homologous nucleotides.

The diagnostic primers for use in the methods of the invention conveniently comprise one or more of

5T variant primer 5'(N)nTGTTAAAGAC3',

7T variant primer 5'(N\*)n\*(N)nTAAAAAAGC3' and

9T variant primer 5'(N\*)n\*(N)nAAAAAAAATC3' wherein N, n, N\* and n\* are as defined above and n is an integer between 6 and 26.

It will be understood that the nucleotide sequence as defined by (N)n in the diagnostic primer is normally selected to be 100% complementary to the corresponding genomic sequence. However, as required, one or more mismatched bases may be included, for example at the 5' terminus of the primer. For example up to two, three, four or five mismatched base pairs may be included in the nucleotide sequence defined by (N)n. It will also be understood that any mismatched bases must not significantly impair the discriminatory properties of the diagnostic primer.

The integer n is for example 10, up to 15, up to 20, up to 25, or up to 30.

Preferred diagnostic primers include

5T variant primer 5'TAATTCCCCAAATCCCTGTAAAGAC3',

7T variant primer 5'(N\*)n\*TAATTCCCCAAATCCCTGTAAAAAAGC3' and

9T variant primer 5'(N\*)n\*TAATTCCCCAAATCCCTGTAAAAAAAATC3' wherein N\* and n\* are as defined above.

Preferred diagnostic 7T and 9T primers are: 7T variant primer 5'GTAAATCATTCAGCTACTACGCACCTAATTCCTGTAAAAAAGC3' and 9T variant primer

5'GACTGTACGATACTCATTATATGAAGTCAGCTACTTACCTATAGAACGCTTGC  
TAGTTTAATTCCCCAAATCCCTGTTAAAAAAAATC3'.

The above primers have been shown to detect their respective alleles reliably and robustly. Each of the primers disclosed above represents a further and independent aspect of the invention.

In many situations, it will be convenient to use a diagnostic primer of the invention with a further amplification primer in one or more cycles of PCR amplification. A convenient example of this aspect is set out in our European patent number EP-B1-0332435. Any convenient amplification primer may be used provided that the resulting amplification products are of a suitable size for separation and analysis. The further amplification primer is conveniently the polyT common primer

GTACATAAAACAAGCATCTATTGAAAATATCTGAC and this is used in combination with the intron 8 5T, 7T and 9T allele-specific primers.

Therefore according to a further aspect of the invention we provide a diagnostic method for the detection of the 5T, 7T and 9T alleles in intron 8 of the human CFTR gene which method comprises contacting a test sample of nucleic acid from an individual with a multiplex of diagnostic primers comprising (i) 5T variant primer 5'(N)nAAAGAC3', (ii) 7T variant primer 5'(N\*)n\*(N)nAAAAGC3' and (iii) 9T variant primer 5'(N\*)n\*(N)nAAAATC3', wherein N represents additional nucleotides which base pair with the corresponding genomic sequence in the respective allele and n is an integer between 10 and 30 and N\* represents additional non-homologous nucleotides which do not base pair with the corresponding genomic sequence in the respective allele and n\* is an integer between 5 and 60, and a common amplification primer in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, and subjecting the mixture to PCR amplification such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a corresponding PCR amplification product.

The above aspects of the invention are referred to as primer set 1 and conveniently illustrated by reference to Table 1 and the disclosure of a specific primer mix 1.

References to primer "mixes" and "sets" are not intended to be limiting and the terms are used throughout the text interchangeably.

It will also be appreciated that we do not want to be limited solely to discrimination between the human CFTR gene intron 8 5T, 7T and 9T alleles by allele associated PCR

product sizes; hence a further definition of N\* is one from a group of labels that can be conjugated to primers whereby an individual label is associated specifically to one allele-specific primer and n\* may be 0 or an integer greater than 1. An example of labels according to this definition of N\* is a group of molecules such as fluorophores for example, (Brown T. and Brown DJS. in Newton CR. (ed): PCR - Essential Data: 1st ed. Wiley; Chichester, pp57-71 (1995)).

The primers may be manufactured using any convenient method of synthesis. Examples of such methods may be found in standard textbooks, for example "*Protocols For Oligonucleotides And Analogues: Synthesis And Properties*," Methods In Molecular Biology Series; Volume 20; Ed. Sudhir Agrawal, Humana ISBN: 0-89603-247-7; 1993; 1<sup>st</sup> Edition.

It will be appreciated that any of the above diagnostic methods may, if appropriate, also be configured so that extension of the diagnostic primer indicates the absence of the respective CFTR gene polyT allele.

As outlined earlier above, the polyT methods of this invention may be used with any known CFTR testing procedure.

In addition we have now devised novel diagnostic primer sequences for the detection of the W1282X, 1717-1G>A, G542X, N1303K, DF508 and 3849+10kbC>T mutations of the human CFTR gene using ARMS allele specific amplification. Therefore the polyT multiplex of the present invention is conveniently accompanied by the use in a separate ARMS reaction of one or more of:

DF508 mutant primer 5'(N)nACCATT3',

3849+10kb C>T mutant primer 5'(N)nTACGCA3',

N1303K mutant primer 5'(N)nTCCATC3',

1717-1G>A mutant primer 5'(N)nTAATTA3',

W1282X mutant primer 5'(N)nCAGTCA3', and

G542X mutant primer 5'(N)nTTCTCT3' wherein N and n are as previously defined, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

Two or more of the above diagnostic primers are conveniently used as a multiplex and more conveniently with a suitable amplification primer. Preferably several or all of the above primer sequences are used in a single multiplex reaction.

Convenient primers include:

W1282X mutant primer 5'(N)nGCAACAGTCA3',

1717-1G>A mutant primer 5'(N)nTTGGTAATTA3',

G542X mutant primer 5'(N)nATAGTTCTCT3',

5 N1303K mutant primer 5'(N)nGGGATCCATC3',

DF508 mutant primer 5'(N)nAAACACCATT3' and

3849+10kb C>T mutant primer 5'(N)nGTCTTACGCA3' wherein N, n, are as defined above and n is an integer between 6 and 26.

Preferred diagnostic primers include

10 W1282X mutant primer 5'TCTTGGGATTCAATAACTTTGCAACAGTCA3',

1717-1G>A mutant primer 5'TCTCGAATTTTCTATTTTGGTAATTA3',

G542X mutant primer 5'AGTTTGCAGAGAAAGACAATATAGTTCTCT3',

N1303K mutant primer 5'TGATCACTCCACTGTTTCATAGGGATCCATC3',

DF508 mutant primer 5'GTATCTATATTCATCATAGGAAACACCATT3', and

15 3849+10kb C>T mutant primer 5'GAACATTTTCCTTTCAGGGTGTCTTACGCA3'.

The above aspects of the invention are referred to as primer set 2A and illustrated by reference to Table 3 and the specific primer mix 2A. Primer mix 2A is conveniently used in combination with primer mix 1 as set out in Table 1.

We have also devised novel diagnostic primer sequences for detection by ARMS allele  
20 specific amplification of the DF508 normal allele, and the 621+1G>T, R117H, R334W, G551D, R553X, and R1162X mutations of the human CFTR gene. Therefore the polyT multiplex of the present invention is conveniently accompanied by the use in a separate ARMS reaction of one or more of:

DF508 non-mutant primer 5'(N)nACCACA3',

25 W1282X mutant primer 5'(N)nCAGTCA3',

1717-1 mutant primer 5'(N)nTAATTA3',

G542X mutant primer 5'(N)nTTCTCT3',

N1303K mutant primer 5'(N)nTCCATC3',

DF508 non-mutant primer 5'(N)nACCACA3',

30 DF508 mutant primer 5'(N)nACCATT3' and

3849+10kb C>T mutant primer 5'(N)nTACGCA3' wherein N and n are as previously defined, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic



variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

Two or more of the above diagnostic primers are conveniently used as a multiplex and more conveniently with a suitable amplification primer. Preferably several or all of the above  
5 primer sequences are used in a single multiplex reaction.

Convenient diagnostic primers include:

DF508 non-mutant primer 5'(N)nAAACACCACA3',

R117H mutant primer 5'(N)nGCGATAGACT3',

621+1G>T mutant primer 5'(N)nGAAGTATTGA3',

10 R334W mutant primer 5'(N)nATCATCCTGT3',

R1162X mutant primer 5'(N)nTCTGTGAGTT3',

R553X mutant primer 5'(N)nTTCTTGCTGA3' and

G551D mutant primer 5'(N)nGCTCGTTGTT3' wherein N, n, are as defined above and n is an integer between 6 and 26.

15 Preferred diagnostic primers include:

R117H mutant primer 5'AGCCTATGCCTAGATAAATCGCGATAGACT3',

621+1G>T mutant primer 5'TGCCATGGGGCCTGTGCAAGGAAGTATTGA3',

R334W mutant primer 5'CCTATGCACTAATCAAAGGAATCATCCTGT3',

R1162X mutant primer 5'TATTTTATTTTCAGATGCGATCTGTGAGTT3',

20 R553X mutant primer 5'TTATTCACCTTGCTAAAGAAATTCTTGCTGA3',

G551D mutant primer 5'GCTAAAGAAATTCTTGCTCGTTGTT3'.

The above aspects of the invention are referred to as primer set 2B and are conveniently illustrated by reference to Table 3 and the specific primer mix 2B. Primer mix 2B is conveniently used in combination with mix1 as set out in Table 1. We have also devised novel  
25 diagnostic primer sequences for the A455E, 2183AA>G, 3659delC, D1507, 1078delT, R347P, S1251N and E60X mutations of the human CFTR gene using ARMS allele specific amplification. Therefore the polyT multiplex of the present invention is conveniently accompanied by the use in a separate ARMS reaction of one or more of:

A455E mutant primer 5'(N)nGTTGTA3',

30 2183AA>G mutant primer 5'(N)nGATAGC3',

3659delC mutant primer 5'(N)nCCTAGA3',

D1507 mutant primer 5'(N)nATAACT3',

- 11 -

1078delT mutant primer 5'(N)nTTCCTG3',  
 R347P mutant primer 5'(N)nTCTACC3',  
 S1251N mutant primer 5'(N)nGAAGCA3' and  
 E60X mutant primer 5'(N)nCAGTTA3'

5 wherein N and n are as previously defined, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

10 Two or more of the above diagnostic primers are conveniently used as a multiplex and more conveniently with a suitable amplification primer. Preferably several or all of the above primer sequences are used in a single multiplex reaction.

The diagnostic primer conveniently comprises one or more of

A455E mutant primer 5'(N)nAGTTGTTGTA3',  
 15 1078delT mutant primer 5'(N)nAGGGTTCCTG3',  
 R347P mutant primer 5'(N)nTTGTTCTACC3',  
 DI507 mutant primer 5'(N)nGAAAATAACT3',  
 3659delC mutant primer 5'(N)nTAAACCTAGA3',  
 2183AA>G mutant primer 5'(N)nAAAAGATAGC3',  
 20 S1251N mutant primer 5'(N)nCAGGGAAGCA3' and  
 E60X mutant primer 5'(N)nAAGCCAGTTA3' wherein N, n, are as defined above and n is an integer between 6 and 26.

Preferred diagnostic primers include

A455E mutant primer 5'TTCAAGATAGAAAGAGGACAGTTGTTGTA3',  
 25 1078delT mutant primer 5'CCTTCTTCTTCTCAGGGTTCCTG3',  
 R347P mutant primer 5'CACCATCTCATTCTGCATTGTTCTACC3',  
 DI507 mutant primer 5'GCCTGGCACCATTAAAGAAAATAACT3',  
 3659delC mutant primer 5'ATGCCAACAGAAGGTAAACCTAGA3',  
 2183AA>G mutant primer 5'CAAACCTCTCCAGTCTGTTTAAAGATAGC3',  
 30 S1251N mutant primer 5'GGAAGAACTGGATCAGGGAAGCA3' and  
 E60X mutant primer 5'TTAGGATTTTCTTTGAAGCCAGTTA3'.

The above aspects of the invention are referred to as primer set 2C and conveniently illustrated by reference to Table 4 and the specific primer mix 2C. Primer mix 2C is conveniently used in combination with mix1 as set out in Table 1.

We have also devised novel diagnostic primer sequences for the detection of the  
5 G85E, 405+1G>A, S549R, W1089X and D1152H mutations of the human CFTR gene using ARMS allele specific amplification. Therefore the polyT multiplex of the present invention is conveniently accompanied by the use in a separate ARMS reaction of one or more of :

G85E mutant primer 5'(N)nCTACGA3',

405+1G>A mutant primer 5'(N)nTAGTGA3',

10 S549R mutant primer 5'(N)nCTGACG3',

W1089X mutant primer 5'(N)nCAAATA3' and

D1152H mutant primer 5'(N)nCACTTG3' wherein N and n are as previously defined, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in  
15 the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

Two or more of the above diagnostic primers are conveniently used as a multiplex and more conveniently with a suitable amplification primer. Preferably several or all of the above primer sequences are used in a single multiplex reaction.

20 The diagnostic primer conveniently comprises one or more of

G85E mutant primer 5'(N)nTGTCTACGA3',

405+1G>A mutant primer 5'(N)nTATTTAGTGA3',

S549R mutant primer 5'(N)nCACTGACG3',

W1089X mutant primer 5'(N)nCTGCCAAATA3',

25 D1152H mutant primer 5'(N)nTATCCACTTG3' wherein N, n, are as defined above and n is an integer between 6 and 26.

Preferred diagnostic primers include:

G85E mutant primer

5'TAGCCATTGATGACGGAGCGATGTTTTCTGGAGATTATGTTCTACGA3'

30 405+1G>A mutant primer

5'GATTTATGTTCTATGGAATCTTTTATATTTAGTGA3',

S549R mutant primer 5'TGGAGAAGGTGGAATCACTGACG3',

W1089X mutant primer 5'AAGCTCTGAATTTACATACTGCCAAATA3' and  
D1152H mutant primer 5'AAAGATGATAAGACTTACCAAGCTATCCACTTG3'

The above aspects of the present invention are referred to as primer set 3 and conveniently illustrated by reference to Table 5 and the specific mix 3. Mix 3 may for example be used in combination with mix1 as set out in Table 1.

The above primers in mixes 2A, 2B, 2C and 3 have been shown to detect their respective alleles reliably and robustly.

The primers may be manufactured using any convenient method of synthesis. Examples of such methods may be found in standard textbooks, for example "*Protocols For Oligonucleotides And Analogues: Synthesis And Properties*," Methods In Molecular Biology Series; Volume 20; Ed. Sudhir Agrawal, Humana ISBN: 0-89603-247-7; 1993; 1<sup>st</sup> Edition.

It will be appreciated that any of the above diagnostic methods may, if appropriate, also be configured so that extension of the diagnostic primer indicates the absence of the respective CFTR gene mutation or polyT allele.

The test sample of nucleic acid is preferably a blood sample but may also conveniently be a sample of any body fluid, or tissue obtained from an individual. The individual is any convenient mammal, preferably a human being. It will be appreciated that the test sample may equally be a nucleic acid sequence corresponding to the sequence in the test sample. That is to say that all or a part of the region in the sample nucleic acid may firstly be amplified using any convenient technique such as PCR before use in the method of the invention.

Any convenient enzyme for polymerisation may be used provided that it does not affect the ability of the DNA polymerase to discriminate between normal and mutant template sequences to any significant extent. Examples of convenient enzymes include thermostable enzymes which have no significant 3'-5' exonuclease activity, for example *Taq* DNA polymerase, particularly "Ampli Taq Gold"<sup>TM</sup> DNA polymerase (PE Applied Biosystems), Stoffel fragment, or other appropriately N-terminal deleted modifications of *Taq* or *Tth* (*Thermus thermophilus*) DNA polymerases.

Significantly the diagnostic methods of the invention can be used in combination to enable the presence or absence of mutations or polyT alleles of the CFTR gene to be detected simultaneously (i.e. the method comprises several multiplex tests). Primer Mix 1 can be used in isolation or in combination with any other cystic fibrosis diagnostic tests. Because the multiplex polyT test using the mix 1 primer mix adds value to the information derived from

other cystic fibrosis diagnostic tests it is convenient to use the mix 1 primer mix in combination with the mix 2A primer mix. Similarly, it is convenient to use the mix 1 primer mix in combination with the mix 2B primer mix. Similarly, it is convenient to use the mix 1 primer mix in combination with the mix 2C primer mix. Similarly, it is convenient to use the mix 1 primer mix in combination with the mix 3 primer mix. Similarly, it is convenient to use the mix 1 primer mix in combination with any two of the mix 2A, 2B, 2C or 3 primer mixes. Similarly, it is convenient to use the mix 1 primer mix in combination with any three of the mix 2A, 2B, 2C or 3 primer mixes. Similarly, it is convenient to use the mix 1 primer mix in combination with all of the mix 2A, 2B, 2C or 3 primer mixes.

The above test combinations give genotype information by distinguishing between individuals who are heterozygous and homozygous for either of the CFTR polyT alleles and the more common DF508 CFTR gene mutation. If desired, the tests may provide genotype information for the 1717-1G>A, G542X, W1282X, N1303K, 3849+10kbC>T, 621+1G>T, R553X, G551D, R117H, R1162X, R334W, A455E, 2183AA>G, 3659delC, 1078delT, DI507, R347P, S1251N, E60X, G85E, 405+1G>A, S549R, W1089X, and D1152H CFTR gene mutations by the inclusion of primers specific for the normal CFTR gene sequences that correspond to the respective mutated CFTR gene sequences associated with these mutations. Whilst each multiplex test may have any of the primers described above in combination the preferred diagnostic primer multiplexes are the mix 1 comprising primers for the 5T, 7T and 9T alleles. Mix 2A comprising primers for the 1717-1G>A, G542X, W1282X, N1303K, DF508 and 3849+10kbC>T mutations; mix 2B comprising primers for the 621+1G>T, R553X, G551D, R117H, R1162X and R334W mutations and the normal CFTR gene sequence corresponding to the DF508 mutation; mix 2C comprising primers for the A455E, 2183AA>G, 3659delC, 1078delT, DI507, R347P, S1251N and E60X mutations; mix 3 comprising primers for the G85E, 405+1G>A, S549R, W1089X, and D1152H mutations.

We have developed a validated test using primer mix 2A for the 1717-1G>A, G542X, W1282X, N1303K, DF508 and 3849+10kbC>T mutations and have applied these in a thorough investigation of the incidence of the mutations in 100 individuals.

In many situations, it will be convenient to use a diagnostic primer of the invention with a further amplification primer in one or more cycles of PCR amplification. A convenient example of this aspect is set out in our European patent number EP-B1-0332435. In the present application, the further amplification primer is either a 1717-1G>A common, G542X common, W1282X common, N1303K common, DF508 common or 3849+10kbC>T common

primer. The 1717-1G>A common primer TAATCTCTACCAAATCTGGATACTATACC is conveniently used in combination with the 1717-1G>A mutant primer, the G542X common primer TAATCTCTACCAAATCTGGATACTATACC is conveniently used in combination with the G542X mutant primer, the W1282X common primer

- 5 GAATTCCCAAACCTTTTAGAGACATC is conveniently used in combination with the W1282X mutant primer, the, N1303K common primer CTTGATGGTAAGTACATGGGTTTTTCTTAT is conveniently used in combination with the N1303K mutant primer, the DF508 common primer CCAGACTTCACTTCTAATTATGATTATGGG is conveniently used in combination with
- 10 the DF508 mutant primer, the 3849+10kbC>T common primer TTGTGGATCAAATTTTCAGTTGACTTGTCATC is conveniently used in combination with the 3849+10kbC>T mutant primer.

- Any convenient control primer may be used. However, the amplification products from the Mix 1 primer mix serve as internal controls and so an additional amplification
- 15 control is not required.

We have developed a validated test using primer mix 2B for the 621+1G>T, R553X, G551D, R117H, R1162X and R334W mutations and the normal DF508 allele and have applied these in a thorough investigation of the incidence of the mutations in 100 individuals.

- In many situations, it will be convenient to use a diagnostic primer of the invention
- 20 with a further amplification primer in one or more cycles of PCR amplification. A convenient example of this aspect is set out in our European patent number EP-B1-0332435. In the present application, the further amplification primer is either a 1717-1G>A common, G542X common, W1282X common, N1303K common, DF508 common or 3849+10kbC>T common primer. The DF508 common primer GACTTCACTTCTAATGATGATTATGGGAGA is
- 25 conveniently used in combination with the DF508 normal primer, the 621+1G>T common primer GTTTCACATAGTGTATGACCCTCTATATACACTCATT is conveniently used in combination with the 621+1G>T mutant primer, the R553X common primer ATCTAAAATTGGAGCAATGTTGTTTTTGACC is conveniently used in combination with the R553X mutant primer, the G551D common primer
- 30 ATCTAAAATTGGAGCAATGTTGTTTTTGACC is conveniently used in combination with the G551D mutant primer, the R117H common primer GTTTCACATAGTGTATGACCCTCTATATACACTCATT is conveniently used in combination with the R117H mutant primer, the R1162X common primer

TTTTGCTGTGAGATCTTTGACAGTCATTT is conveniently used in combination with the R1162X mutant primer, the R334W common primer

TTTGTATTGCTCCAAGAGAGTCATACCA is conveniently used in combination with the R334W mutant primer.

- 5        Any convenient control primer may be used. We have selected control primers from unrelated regions of the genome, namely, part of the human apolipoprotein B gene and part of the ornithine decarboxylase gene.

- We have developed a validated test using primer mix 2C for the A455E, 2183AA>G, 3659delC, 1078delT, DI507, R347P, S1251N and E60X mutations and have applied these in a  
10 thorough investigation of the incidence of the mutations in 100 individuals.

- In many situations, it will be convenient to use a diagnostic primer of the invention with a further amplification primer in one or more cycles of PCR amplification. A convenient example of this aspect is set out in our European patent number EP-B1-0332435. In the present application, the further amplification primer is either a 621+1G>T common, R553X  
15 common, G551D common, R117H common, R1162X common, R334W common or A455E common primer. The A455E common primer

- GACTGACTGACTGACTGAAATGGAGACTTTTTGTTTATGTGGTTACTAA is conveniently used in combination with the A455E mutant primer, the 2183AA>G common primer GTATGATAGAGATTATATGCAATAAAACATTAACA is conveniently used in  
20 combination with the 2183AA>G mutant primer, the 3659delC common primer TGTGTCTAATATTGATTCTACTGTACAATAATAA is conveniently used in combination with the 3659delC mutant primer, the 1078delT common primer

- ATTTTCCAACTTCATTAGAACTGATCTATTGAC is conveniently used in combination with the 1078delT mutant primer, the DI507 common primer  
25 CACAGTAGCTTACCCATAGAGGAAACA is conveniently used in combination with the DI507 mutant primer, the R347P common primer

- ATTTTCCAACTTCATTAGAACTGATCTATTGAC is conveniently used in combination with the R347P mutant primer, the S1251N common primer  
GCTCACCTGTGGTATCACTCCAA is conveniently used in combination with the S1251N  
30 mutant primer, the E60X common primer

AATCAAACCTATGTTAAGGGAAATAGGACAACTAA is conveniently used in combination with the E60X mutant primer.

Any convenient control primer may be used. We have selected control primers from unrelated regions of the genome, namely, part of the human apolipoprotein B gene and part of the ornithine decarboxylase gene.

We have developed a validated test using primer mix 3 for the G85E, 405+1G>A, S549R, W1089X, and D1152H mutations and have applied these in a thorough investigation of the incidence of the mutations in a number of individuals.

In many situations, it will be convenient to use a diagnostic primer of the invention with a further amplification primer in one or more cycles of PCR amplification. A convenient example of this aspect is set out in our European patent number EP-B1-0332435. In the present application, the further amplification primer is either a G85E common, 405+1G>A common, S549R common, W1089X common, D1152H common primer. The G85E and the 405+1G>A common primer CGATTCGATTCAGTTTCTGTGGTTTCTTAGTGTTTGGA is conveniently used in combination with the G85E mutant primer and/or the 405+1G>A mutant primer, the S549R common primer GTAATTTTTTTACATGAATGACATTTACAGCAA is conveniently used in combination with the S549R mutant primer, the W1089X common primer GGAAATTATTTGTTTAACAATAAAACAATGGAA is conveniently used in combination with the W1089X mutant primer and the D1152H common primer CCAACAACACCTCCAATACCAGTAAC is conveniently used in combination with the D1152H mutant primer.

Any convenient control primer may be used. We have selected control primers from unrelated regions of the genome, namely, part of the human apolipoprotein B gene and part of the ornithine decarboxylase gene.

It will be appreciated that by combining intron 8 polyT allele-specific primers and their respective common primer and/or different mutant or normal primers and their respective common primers, or by using combinations of the 1 and 2A and/or 2B and/or 2C and/or 3 primer sets or mixes for a given DNA sample, further useful diagnostic tests are provided which permit the simultaneous detection of alternative polyT alleles and/or several mutations in the CFTR gene. Each of these combinations will also include the addition of the appropriate control primers as set out in Tables 1 to 5 below. The combination of different primers for the simultaneous detection of two or more point mutations is termed "multiplexing" (see EP-B1-0332435).



A variety of methods may be used to detect the presence or absence of diagnostic primer extension products and/or amplification products. These will be apparent to the person skilled in the art of nucleic acid detection procedures. Preferred methods avoid the need for radiolabelled reagents. Particular detection methods include size separation of amplification products, for example as described in our patents numbers EP 0 332 435 & US 5595890, "ALEX" product detection (Haque K. et al. *Diag. Mol. Pathol.* &: 248-252 (1998)), the detection of amplification incorporated "Sunrise" probes (Nazarenko IA. et al. *Nucl. Acids Res.* 25: 2516-2521 (1997)), the detection of amplification incorporated "Scorpions" primers (Whitcombe D. et al. *Nature Biotechnol.* 17: 804-807 (1999)), the detection of fluorescence polarisation signal combined with dual labelled ARMS primers (Gibson NJ. et al. *Clin. Chem.* 43: 1336-1341 (1997)) patent number EP 0 382 433, the detection of intercalating dyes into ARMS products (Brownie J. et al. *Nucl. Acids Res.* 25: 3235-3241 (1997)) patent number EP 0 731 177, "taqman" product detection for example as described in patent numbers US-A-5487972 & US-A-5210015; and "Molecular Beacons" product detection outlined in patent number WO-95/13399.

One or more of the primer mixes of the invention may be conveniently packaged with instructions for use in the method of the invention and appropriate packaging and sold as a kit. Convenient primer mixes include intron 8 polyT allele-specific primers and CFTR mutation-specific primers, equivalent normal-specific primers; all as hereinbefore disclosed. The kits will conveniently include one or more of the following: appropriate nucleotide triphosphates, for example one or more of dATP, dCTP, dGTP, and dTTP, a suitable polymerase as previously described, and a convenient buffer solution.

The invention will now be illustrated but not limited by reference to the following Examples, Tables and Figures in which:

Table 1 shows a specific primer mix 1

Table 2 shows a specific primer mix 2A

Table 3 shows a specific primer mix 2B

Table 4 shows a specific primer mix 2C

Table 5 shows a specific primer mix 3

Figure 1 shows the relative location of the PCR product bands on an agarose electrophoresis gel corresponding to (PolyT)

Figure 2 shows the relative location of the PCR product bands on an agarose electrophoresis gel corresponding to a CFTR DF508 heterozygote, DF508 homozygote, DF508/ DI507 compound heterozygote, G542X heterozygote, G551D heterozygote, 1078delT heterozygote. (i.e figure from brochure)

- 5        Figure 3 shows the relative location of the PCR product bands on an agarose electrophoresis gel corresponding to a (Middle Eastern Panel)

Figure 4 shows diagrammatically the relative sizes in base pairs and the relative location of the PCR product bands on an agarose electrophoresis gel corresponding to the CFTR intron 8 5T, 7T and 9T alleles.

- 10       Figure 5 shows diagrammatically the size in base pairs and the relative location of the PCR product bands on an agarose electrophoresis gel corresponding to the CFTR 1717-1G>A, G542X, W1282X, N1303K, DF508, 3849+10kbC>T, 621+1G>T, R553X, G551D, R117H, R1162X, R334W, A455E, 2183AA>G, 3659delC, 1078delT, DI507, R347P, S1251N, E6OX mutant alleles and the non-DF508 allele.

- 15       Figure 6 shows diagrammatically the size in base pairs and the relative location of the PCR product bands on an agarose electrophoresis gel corresponding to the CFTR D1152H, W1089X, G85E, 405+1G>A, S549R mutant alleles.

**Examples****Example 1.1 Materials Provided in a diagnostic kit containing 50tests (polyT)**

1. 50 Vial 1 (colour coded) containing primers (5T, 7T, 9T), and deoxynucleotide  
5 triphosphates in buffer.
2. 1 vial dilution buffer.
3. 1 vial AmpliTaq Gold.
4. 1 vial loading dye.
5. 1 vial DNA control, contains human DNA in buffer.
- 10 6. Instructions for use.

**Example 2.1 Materials Provided in a diagnostic kit containing 50 tests (CF20)**

1. 50 Vial 2A (colour coded) containing primers (mutant 1717-1G>A, G542X,  
W1282X, N1303K, DF508 and 3849+10kbC>T), control primers and deoxynucleotide  
15 triphosphates in buffer.
2. 50 Vial 2B (colour coded) containing primers (mutant 621+1G>T, R553X,  
G551D, R117H, R1162X, R334W, normal DF508), control primers and deoxynucleotide  
triphosphates in buffer.
3. 50 Vial 2C (colour coded) containing primers (mutant A455E, 2183AA>G,  
20 3659deIC, 1078deIT, DI507, R347P, S1251N, E6OX), control primers and deoxynucleotide  
triphosphates in buffer.
4. 1 vial dilution buffer.
5. 1 vial AmpliTaq Gold.
6. 1 vial loading dye.
- 25 7. 1 vial normal DNA control, contains human DNA unaffected by the  
mutations detected by the kit, in buffer.
8. Instructions for use.

**Example 3.1 Materials Provided in a diagnostic kit containing 50tests (Middle Eastern)**

- 30 1. 50 Vial 3 (colour coded) containing primers (D1152H, W1089X, G85E,  
405+1G>A, S549R ), control primers and deoxynucleotide triphosphates in buffer
2. 1 vial dilution buffer.

3. 1 vial AmpliTaq Gold.
4. 1 vial loading dye.
5. 1 vial normal DNA control, contains human DNA unaffected by the mutations detected by the kit, in buffer.
- 5 6. Instructions for use.

**Example 4.1 Instructions for CF20 PCR amplification procedure**

1. Program the thermal cycler for an activation program which holds the vials at 94 °C for 20 minutes and an amplification program of 0.5 minutes at 94 °C (denaturation), 2 minutes at 58 °C (annealing) and 1 minute at 72 °C (extension) for 35 cycles. Followed by a 20 minute 72 °C extension stage..
2. Label one vial A, one vial B and one vial C for each sample and control.
3. Microfuge vials A, B and C until all liquid is at the bottom of each vial.
4. Prepare sufficient dilution of the AmpliTaq Gold for the number of samples and controls to be tested. For 10 samples or controls pipette 68 µL sterile deionized water, 20 µL dilution buffer, 100 µL loading dye and 12 µL AmpliTaq Gold into a microfuge tube and mix gently.
5. Carefully open the vial lid and pipette 5 µL of the enzyme dilution into the A, B and C vials using separate tips and re-cap.
- 20 6. Pipette 5 µL of test or Normal DNA Control sample to each of a vial A, B and C vials using separate tips. Add one drop of Sigma light white mineral oil to cover the aqueous phase if using a thermal cycler without a heated lid. Re-cap firmly.
7. For the negative control add 5 µL sterile water to a vial of each of A, B and C. Add one drop of Sigma light white mineral oil to cover the aqueous phase if using a thermal cycler without a heated lid. Re-cap firmly.
- 25 8. Microfuge vials A, B and C until all liquid is at the bottom of each vial.
9. Place all tubes firmly in the thermal cycler block. Initiate the 94 °C activation program. On completion of the activation program, run the amplification program.
10. On completion of the final extension stage, the samples may be stored at room temperature overnight or at 2-8 °C for up to 7 days before analysis by gel electrophoresis.
- 30

**Instructions for CF Poly-T and CF-MEP PCR amplification procedure**

1. Program the thermal cycler for an activation program which holds the vials at 94 °C for 20 minutes and an amplification program of 0.5 minutes at 94 °C (denaturation), 2 minutes at 58 °C (annealing) and 1 minute at 72 °C (extension) for 35 cycles. Followed by a 20 minute 72 °C extension stage..

5 2. Label a vial for each sample and control.

3. Microfuge vials until all liquid is at the bottom of each vial.

4. Prepare sufficient dilution of the AmpliTaq Gold for the number of samples and controls to be tested. For 10 samples or controls pipette 68 µL sterile deionized water, 20 µL dilution buffer, 100 µL loading dye and 12 µL AmpliTaq Gold into a microfuge tube and  
10 mix gently.

5. Carefully open the vial lid and pipette 5 µL of the enzyme dilution into the vials using separate tips and re-cap.

6. Pipette 5 µL of test or DNA Control sample to each vial using separate tips. Add one drop of Sigma light white mineral oil to cover the aqueous phase if using a thermal  
15 cycler without a heated lid. Re-cap firmly.

7. For the negative control add 5 µL sterile water to a separate vial. Add one drop of Sigma light white mineral oil to cover the aqueous phase if using a thermal cycler without a heated lid. Re-cap firmly.

8. Microfuge vials until all liquid is at the bottom of each vial.

20 9. Place all tubes firmly in the thermal cycler block. Initiate the 94 °C activation program. On completion of the activation program, run the amplification program.

10. On completion of the final extension stage, the samples may be stored at room temperature overnight or at 2-8 °C for up to 7 days before analysis by gel electrophoresis.

#### 25 Example 4.2 Procedure

1. 15 x 12 cm horizontal submarine gels with combs of 1.5 mm x 5 mm suspended 1 mm above the gel tray, were prepared using 100 mL of 3% NuSieve™ (FMC Corporation) 3:1 agarose in 134 mM (16.2 g/L) Tris-base, 74.9 mM (4.63 g/L) boric acid, 2.55 mM (0.95 g/L) EDTA buffer with 0.1 µg/mL ethidium bromide (TBE/EtBr). TBE/EtBr was  
30 also used as the running buffer.

2. A 50 Base-Pair Ladder (Amersham Pharmacia Biotech) at 1.5 µg/15µL was prepared in the loading dye (80 µL distilled water / 10µL loading dye / 10µL 50 Base-Pair Ladder). 15µL of this dilution was run adjacent to samples as a molecular weight marker.

3. 15 µL of PCR products from samples or control vials were loaded on a gel (for CF20 A, B and C trios were loaded in adjacent wells).

4. Electrophoresis was carried out at 5 to 6 V/cm between electrodes until the dye front had migrated 5 cm from the loading wells towards the anode (1.5 to 2 hours).

5. After electrophoresis the gels were placed on a UV transilluminator at 260 nm then visualised and photographed.

#### Example 4.3 Interpretation of Results for CF20 and CF-MEP

PCR products will be observed as bands in the vial tracks of the gel.

1. The upper and lower control band must be clearly visible in all samples.  
2. All tracks should be free of excessive smearing and background fluorescence.  
3. The position of the upper and lower control bands should indicate the correct molecular size (see Figure 5 or Figure 6).

4. The negative control should show no bands in the tracks within the area corresponding to the upper and lower control bands. A diagnostic band should not be interpreted if a similar band is also seen in the negative control for that PCR run as this is indicative of contamination with genomic DNA.

If any of the above points are not observed the results should not be interpreted and a repeat test carried out.

5. An individual has two copies of the CFTR gene. Where these copies have the same sequence for any given site, an individual is described as being homozygous for this site. Where the copies differ in sequence at a given site, an individual is described as being heterozygous for this site.

6. The presence of PCR product generated from the normal DF508 primer in vial B of CF20 indicates that the sample contains the normal sequence for this site. The normal PCR product will be observed in the vial B track of the gel at 160bp and is identified by comparison of the band position with an adjacent marker track.

7. For CF20 the PCR products from an individual carrying any of the 1717-IG>A, G542X, W1282X, N1303K, DF508 and 3849+10kbC>T mutations will be observed in the

vial A track of the gel and are identified by comparison of the band position with an adjacent marker track. The product band sizes in base pairs are shown in Figure 5. Only product bands of the correct size should be interpreted.

8. For CF20 the PCR products from an individual carrying any of the 621+1G>T, R553X, G551D, R117H, R1162X, and R334W mutations will be observed in the vial B track of the gel and are identified by comparison of the band position with an adjacent marker track. The product band sizes in base pairs are shown in Figure 5. Only product bands of the correct size should be interpreted.

9. For CF20 the PCR products from an individual carrying any of the A455E, 2183AA>G, 3659delC, 1078delT, D1507, R347P, S1251N and E60X mutations will be observed in the vial C track of the gel and are identified by comparison of the band position with an adjacent marker track. The product band sizes in base pairs are shown in Figure 5. Only product bands of the correct size should be interpreted.

#### 15 Interpretation of Results for CF Poly-T

PCR products will be observed as bands in the vial tracks of the gel.

1. The bands should be of similar intensity to the 250bp band in the 50bp ladder, loaded at 1.5mg/15ml (Amersham Pharmacia Biotech).

2. At least one of the CF Poly-T products should be visible in each DNA sample.

20 3. No bands should be visible in the negative control

4. If any of the above points are not observed the results should not be interpreted and a repeat test carried out.

5. For CF Poly-T the product band sizes in base pairs are shown in Figure 4. Only product bands of the correct size should be interpreted.

25 6. As a result of allele expansion within the microsatellite region (TG)<sub>8-13</sub>, contiguous to the Poly-T region, the length of the diagnostic PCR products may vary. As a result very small product band shifts may be observed during electrophoresis. These shifts do not alter the interpretation of test results.

7. Owing to the nature of the sequence surrounding the Poly-T locus, 30 heteroduplexes may occasionally be formed from the PCR products and may therefore be visible in some test sample tracks. The position of these heteroduplexes is indicated in figure

8. The observed heteroduplexes are summarised as follows:

A 7T/7T individual may have a heteroduplex band larger than the 7T allele but not as large as the 9T allele.

A 5T/7T individual may have a heteroduplex that runs in a similar position to the 9T allele. Therefore if all three diagnostic bands are visible the individual's genotype is 5T/7T.

- 5 The heteroduplex, if present, is generally weaker than the corresponding 5T and 7T diagnostic products.

A 5T/9T individual occasionally has a heteroduplex which run larger than the 9T allele.

- 10 Owing to the sequence variability in the region proximal to the Poly-T repeat, additional (previously unseen) heteroduplexes may be observed.

#### Example 4.4 Performance Characteristics

- For CF20 one hundred EDTA blood samples were tested using the procedures described in Examples 4.1 and 4.2 in a 'blind' in-house study. Samples were prepared using the method described in Example 4.5. Each result obtained was confirmed for each of the 20 mutations using alternative methods. Of the 100 individuals tested, 92 were normal and 8 were DF508 heterozygotes.

- 20 In addition blood and mouthwash samples were taken from each of 40 individuals and tested using the diagnostic assay. The result obtained from each mouthwash sample was concordant with that obtained using the blood sample from the same individual.

A number of compound heterozygotes within the same multiplex test were sourced and analysed using the diagnostic assay. In all cases, both mutations were detected. Those tested were:

The 'A mix'

- 25 1717-1G>A / N1303K, 1717-1G>A / DF508, G542X / W1282X, G542X / N1303K, G542X / DF508, W1282X / N1303K, W1282X / DF508, W1282X / 3849+10kbC>T, N1303K / DF508, DF508 / 3849+10kbC>T.

The 'B mix'

R553X / G551D, R553X / R334W.

- 30 Owing to the rarity of the remaining combinations and the subsequent unavailability of samples, other compound heterozygotes were not evaluated.



The rare R1283M, 1717-2A>G, R117C, 621+2T>C, R117L, I506V mutations have been evaluated for cross reactivity. These were tested using the procedures described in Examples 1.2 and 1.3 and were not detected. In addition the following polymorphisms were not detected by the test: 3617G/T, 1655T/G (F508C), 1651A/G.

5 Evaluation of known mutations and polymorphisms in the CFTR gene using the procedures described in Examples 4.1 and 4.2 has highlighted the following observations:

1. DI507 in combination with dupl 716+51 >61 will produce a PCR product of 233bp in the 'C mix' rather than the expected 222bp.
2. Slight cross reactivity of the R347P primer in the 'C mix' with the rare R347H  
10 mutation was observed which results in faint PCR product visible at the R347P position in the 'C mix'.
3. A recently reported mutation, 2184insG, will theoretically cross react with the 2183AA>G primer in the 'C mix'. A diagnostic 2183AA>G band may indicate that the 2184insG mutation is present.
- 15 4. The following mutations, which have not been checked due to unavailability of relevant samples, may interfere with test function: R117P, 621+2T>G, R553G, R553Q, R347L, I506T, I506S and the rare combination of DI507 with the polymorphism 1651A/G. For CF Poly-T one hundred EDTA blood samples were tested using the procedures described in Examples 4.1 and 4.2 in a 'blind' in-house study. Samples were prepared using the method  
20 described in Example 4.5. The genotype for each sample was independently confirmed. Additionally 100 mouthwash samples paired with the blood samples were tested on the CF Poly-T test. The result obtained from each mouthwash sample was concordant with that obtained using the blood sample from the same individual.

There are a number of sequence variants in the CFTR gene of unknown frequency  
25 which are located close to the Poly-T tract. The following mutations, which have not been checked due to unavailability of relevant samples, may theoretically interfere with test function:

1. The 1342-11TTT>G mutation which results in a (TG)<sub>13</sub>(T)<sub>3</sub> individual with an apparent 5T disease state will produce diagnostic PCR product from the 5T primer.  
30
2. The 1342-13G/T polymorphism which alters the final TG repeat to TT by substitution of the final does not affect the performance of the test as this polymorphism changes the status of individuals with respect to the polythymidine tract length.

3. The effects of the 1342-1G>C, and 1342-2A>C and 1342-2delAG mutations have not been tested but may affect the performance of the ELUCIGENE™ CF Poly-T test.

For CF-MEP one hundred EDTA blood samples were tested using the procedures described in Examples 4.1 and 4.2 in a 'blind' in-house study. Samples were prepared using the method described in Example 4.5. All 100 samples tested were normal with respect to the CF-MEP mutations.

One compound heterozygotes within the same multiplex test was sourced and analysed using the diagnostic assay. In the sample G85E/D1152H both mutations were detected

#### 10 Example 4.5 Method for Preparation of DNA from Whole Blood (EDTA) Samples:

1. Pipette 80 µL of each blood sample into a screw-topped microfuge tube.
2. Pipette 320 µL of 170 mM (9.09 g/L) NH<sub>4</sub>Cl solution into each tube.
3. Mix for 20 minutes by gentle swirling and inversion. Avoid vigorous agitation and formation of foam.
4. Centrifuge each tube for 2 minutes at 12 000g until a cell pellet is formed.
5. Using a pipette remove and discard the supernatant liquid.
6. Pipette 300 µL of 10 mM (0.58 g/L) NaCl/10 mM (3.72 g/L) EDTA into each tube and resuspend the cells by vortex mixing.
7. Centrifuge each tube for 1 minute at 12 000g until a cell pellet is formed.
8. Repeat steps 5 to 7 at least a further two times until all visible red coloration in the supernatant liquid has been removed.
9. Using a pipette remove and discard the supernatant liquid.
10. Pipette 200 µL of 50 mM (2 g/L) NaOH solution into each tube and resuspend the cells by vortex mixing.
11. Incubate at 100°C for 10 minutes.
12. Pipette 40 µL of 1 M (121.1 g/L) Tris/HCl (pH 7.5) into each tube and vortex mix.
13. Add 1ml sterile deionised water to each microfuge tube to give a total DNA sample volume of 1.24mL.
14. Centrifuge each tube for 1 minute at 12 000g until a pellet of cell debris is formed. The DNA is contained within the supernatant liquid.

**Example 4.6 Method for Preparation of DNA from Mouthwash Samples:**

1. Agitate 10 mL of 0.9% saline in the mouth for 20 seconds. Collect the suspension in a sterile plastic universal tube.
- 5 2. Pellet the cells by centrifugation at 800g for 10 minutes at 18-28 °C.
3. Using a pipette remove and discard the supernatant liquid.
4. Pipette 500 µL of 10 mM (0.58 g/L) NaCl/10 mM (3.72 g/L) EDTA into each tube and resuspend the cells by vortex mixing.
5. Transfer each sample to a screw-capped microfuge tube.
- 10 6. Centrifuge each tube for 1 minute at 12 000g until a cell pellet is formed.
7. Using a pipette carefully remove and discard the supernatant liquid.
8. Pipette 500 µL of 50 mM (2g/L) NaOH solution into each tube and resuspend the cells by vortex mixing.
9. Incubate at 100°C for 10 minutes
- 15 10. Pipette 100 µL of 1 M (121.1 g/L) Tris/HCl (pH 7.5) into each tube and vortex mix.
11. Centrifuge each tube for 1 minute at 12 000g until a pellet of cell debris is formed. The DNA is contained within the supernatant liquid.
12. Transfer 100 µL of the supernatant (DNA sample) to a fresh, labelled
- 20 microfuge tube.
13. Add 400 µL of sterile deionised water to each DNA sample to give a total volume of 500 µL.

The mix formulation for mix 1 is presented below (Table 1). The 25µl total reaction volume also includes 1x ARMS, 100mM dNTP, 1.5U AmpliTaq Gold, 1x CRS.

WO 01/21833

PCT/GB00/03597

- 29 -

Table 1

Primer	F/R	Primer Ref	Mismatch	Length	Sequence
5T	R	5TMR	GT-3	26	TAA TTC CCC AAA TCC CTG TTA AAG AC
7T	R	7T25	GT-2	53	GTT AAT CAT TCA GCT ACT ACG CAC CTA ATT CCC CAA ATC CCT GTT AAA AAA GC
9T	R	9T60	TT-2	90	GAC TGT ACG ATA CTC ATT TAT ATG AAG TCA GCT TAC TTA CCT ATA GAA CGC TTG CTA GTT TAA TTC CCC AAA TCC CTG TTA AAA AAA ATC
Intron 8 polyT C	F	PT COM FOR		35	GTA CAT AAA ACA AGC ATC TAT TGA AAA TAT CTG AC

The mix formulation for the 2A mix is presented below (Table 2). The 25µl total reaction volume also includes 1x ARMS, 100mM dNTP, 1.5U AmpliTaq Gold, 1x CRS

Table 2

5

Primer	F/R	Primer Ref	Mismatch	Length	Conc µM	Sequence
Apo B	F	ABF	CA-2	23	0.2	GAG CAC AGT ACG AAA AAC CAC CT
Apo B	R	ABR20 A	TG-2 AG-6	25	0.2	AAA CAC GAA GAT GCT GTC TAC TAT C
ODC (FPLC purified)	F	ODFFP LC	CT-2	30	0.3	AGA GGA TTA TCT ATG CAA ATC CTT GTA ACC
ODC (FPLC purified)	R	ODRFP LC	AC-2	30	0.3	TCA ACT TCA CTA TCA AAA GTC ATC ATC TAG
W1282X M	F	WXMF	CC-2	30	1	TCT TGG GAT TCA ATA ACT TTG CAA CAG TCA
W1282X C	R	WXCR		25	1	GAA TTC CCA AAC TTT TAG AGA CAT C
1717-1G>A M	F	17MF20	TT-2 GA-23	40	2	TAC TAA AAG TGA CTC TCG AAT TTT CTA TTT TTG GTA ATT A
Exon 11 C	R	Aex11C R		29	1	TAA TCT CTA CCA AAT CTG GAT ACT ATA CC
G542X M	F	GXMF	CA-2	30	1	AGT TTG CAG AGA AAG ACA ATA TAG TTC TCT

WO 01/21833

PCT/GB00/03597

- 31 -

N1303K C	F	NKCF2 0	TC-9	30	1	CTT GAT GGT AAG TAC ATG GGT TTT TCT TAT
N1303K M	R	NKMR	TT-2	30	1	TGA TCA CTC CAC TGT TCA TAG GGA TCC ATC
DF508 C	F	ADFCF 20	TC-12	30	0.5	CCA GAC TTC ACT TCT AAT TAT GAT TAT GGG
DF508 M	R	ADFMR	TT-2	30	0.5	GTA TCT ATA TTC ATC ATA GGA AAC ACC ATT
3849+10kb C	F	38CF		31	0.2	TTG TGG ATC AAA TTT CAG TTG ACT TGT CAT C
3849+10kb M	R	38MR20	GA-3	29	0.2	GAA CAT TTC CTT TCA GGG TGT CTT ACG CA

The mix formulation for the 2B mix is presented below (Table 3). The 25µl total reaction volume also includes 1x ARMS, 100mM dNTP, 1.5U AmpliTaq Gold, 1x CRS.

5

Table 3

Primer	F/R	Primer Ref	Mismatch	Length	Conc µM	Sequence
Apo B	F	ABF	CA-2	23	0.1	GAG CAC AGT ACG AAA AAC CAC CT
Apo B	R	ABR20 B	CA-2	26	0.1	CAT TTA GTT TCA GCC CAG GAA TAA CG

WO 01/21833

PCT/GB00/03597

- 32 -

ODC (FPLC purified)	F	ODFFP LC	CT-2	30	0.3	AGA GGA TTA TCT ATG CAA ATC CTT GTA ACC
ODC (FPLC purified)	R	ODRFP LC	AC-2	30	0.3	TCA ACT TCA CTA TCA AAA GTC ATC ATC TAG
621+1G>T M	R	62MR	GT-2	30	2	TGC CAT GGG GCC TGT GCA AGG AAG TAT TGA
R117H M	R	RHMR	CC-2	30	0.5	AGC CTA TGC CTA GAT AAA TCG CGA TAG ACT
621/R117H C	F	ex4CF	CT-8 TC- 26 GA-27	37	1.5	GTT TCA CAT AGT GTA TGA CCC TCT ATA TAC ACT CAT T
R334W M	F	RWMF	GG-2	30	0.2	CCT ATG CAC TAA TCA AAG GAA TCA TCC TGT
R334W C	R	RWCR		30	0.2	TTT GTT TAT TGC TCC AAG AGA GTC ATA CCA
G551D M	R	GDMR2 0	TT-2	27	2	TTG CTA AAG AAA TTC TTG CTC GTT GTT
R553X M	R	R5MR	GG-2	60	2	GAC TGA CTG ACT GAC TGA CTC TGA CTG ACT TAT TCA CCT TGC TAA AGA AAT TCT TGC TGA
G551D/ R553X C	F	Bex11C F20	GG-20 GA- 21	31	2	ATC TAA AAT TGG AGC AAT GTT GTT

WO 01/21833

PCT/GB00/03597

- 33 -

						TTT GAC C
DF508 N	R	BDFNR	CT-2	30	0.5	GTA TCT ATA TTC ATC ATA GGA AAC ACC ACA
DF508 C	F	BDFCF		30	0.5	GAC TTC ACT TCT AAT GAT GAT TAT GGG AGA
R1162X M	F	R1MF	TG-2	30	0.2	TAT TTT TAT TTC AGA TGC GAT CTG TGA GTT
R1162X C	R	R1CR		29	0.2	TTT TGC TGT GAG ATC TTT GAC AGT CAT TT

The mix formulation for the 2C mix is presented below (Table 4). The 25µl total reaction volume also includes 1x ARMS, 100mM dNTP, 1.5U AmpliTaq Gold, 1x CRS.

5

Table 4

Primer	F/R	Primer Ref	Mismatch	Length	Conc µM	Sequence
Apo B	F	ABF	CA-2	23	0.2	GAG CAC AGT ACG AAA AAC CAC CT
Apo B	R	ABR20C		30	0.2	CAG CTT CTT ATA GAT TTG TAT TTC TCT GAA
ODC (FPLC purified)	F	ODFFPL C	CT-2	30	0.3	AGA GGA TTA TCT ATG CAA ATC CTT GTA ACC
ODC (FPLC	R	ODRCFP LC	AC-2	31	0.3	TCA ACT TCA CTA TCA AAA GTC ATC



WO 01/21833

PCT/GB00/03597

- 34 -

purified)						ATC TGA A
A455E M	F	A45MF	TC-2	50	1	GAC TGA CTG ACT GAC TGA AAT TTC AAG ATA GAA AGA GGA CAG TTG TTG TA
A455E C	R	A45CR		49	1	GAC TGA CTG ACT GAC TGA AAT GGA GAC TTT TTG TTT ATG TGG TTA CTA A
2183AA>G C	F	2183CF		35	2	GTA TGA TAG AGA TTA TAT GCA ATA AAA CAT TAA CA
2183AA>G M	R	2183MR	AA-3	31	2	CCC AAA CTC TCC AGT CTG TTT AAA AGA TAG C
3659delC M	F	3659MF	GG-2	27	2	GAC ATG CCA ACA GAA GGT AAA CCT AGA
3659delC C	R	3659CR	CA-29	34	2	TGT GTC TAA TAT TGA TTC TAC TGT ACA ATA ATA A
DI507 M	F	DIMF	AA-3	26	0.5	GCC TGG CAC CAT TAA AGA AAA TAA CT
DI507 C	R	DICR		27	0.5	CAC AGT AGC TTA CCC ATA GAG GAA ACA
1078delT M	F	1078MF	CA-3	23	2	CCT TCT TCT TCT CAG GGT TCC TG

WO 01/21833

PCT/GB00/03597

- 35 -

R347P M	F	R34PMF	AC-3	27	0.2	CAC CAT CTC ATT CTG CAT TGT TCT ACC
Exon 7 C	R	ex7CR	CC-24 CC-29	35	1	ATT TTT CCA AAC TTC ATT AGA ACT GAT CTA TTG AC
S1251N M	F	S12MF	CT-2	23	0.2	GGA AGA ACT GGA TCA GGG AAG CA
S1251N C	R	S12CR		23	0.2	GCT CAC CTG TGG TAT CAC TCC AA
E60X C	F	E60XCF		34	0.5	AAT CAA ACT ATG TTA AGG GAA ATA GGA CAA CTA A
E60X M	R	E60XMR	TG-3	26	0.5	TTA GGA TTT TTC TTT GAA GCC AGT TA

The mix formulation for mix 3 is presented below (Table 5). The 25µl total reaction volume also includes 1x ARMS, 100mM dNTP, 1.5U AmpliTaq Gold, 1x CRS.

5

Table 5

Primer	F/R	Primer Ref	Mismatch	Length	Conc µM	Sequence
Apo B	F	ABF	CA-2	23	0.1	GAG CAC AGT ACG AAA AAC CAC CT
Apo B	R	ABR20B	CA-2	26	0.1	CAT TTA GTT TCA GCC CAG GAA TAA CG

WO 01/21833

PCT/GB00/03597

- 36 -

ODC (FPLC purified)	F	ODFFPLC	CT-2	30	0.3	AGA GGA TTA TCT ATG CAA ATC CTT GTA ACC
ODC (FPLC purified)	R	ODRCFPL C	AC-2	31	0.3	TCA ACT TCA CTA TCA AAA GTC ATC ATC TGA A
S549R C	R	122-99		33	0.5	GTA ATT TTT TTA CAT GAA TGA CAT TTA CAG CAA
S549R T>G M	F	126-99	CC-2	25	0.5	TGG AGA AGG TGG AAT CAC ACT GAC G
W1089X C	R	F10682		33	1.0	GGA AAT TAT TTG TTT AAC AAT AAA ACA ATG GAA
W1089X M	F	140-99	AG-3	28	1.0	AAG CTC TGA ATT TAC ATA CTG CCA AAT A
D1152H C	F	F10683		26	0.5	CCA ACA ACA CCT CCA ATA CCA GTA AC
D1152H M	R	145-99	TT-3	33	0.5	AAA GAT GAT AAG ACT TAC CAA GCT ATC CAC TTG
exon3 C (NH-tail)	R	670-99	TC-21, CT-22	39	0.5	CGA TTC GAT TCA GTT TTC TGT GGT TTC TTA

						GTG TTT GGA
G85E M (NH-tail)	F	616-99	CA-3	50	0.5	TAG CCA TTG ATG ACG GAG CGA TGT TTT TTC TGG AGA TTT ATG TTC TAC GA
405+1G>A M	F	134-99	TC-3	36	1.0	GAT TTA TGT TCT ATG GAA TCT TTT TAT ATT TAG TGA

**Claims:**

1. A diagnostic method for the detection of the 5T, 7T and 9T alleles in intron 8 of the human CFTR gene which method comprises contacting a test sample of nucleic acid from an individual with a multiplex of diagnostic primers comprising (i) 5T variant primer 5'(N)nAAAGAC3', (ii) 7T variant primer 5'(N\*)n\*(N)nAAAAGC3' and (iii) 9T variant primer 5'(N\*)n\*(N)nAAAATC3', wherein N represents additional nucleotides which base pair with the corresponding genomic sequence in the respective allele and n is an integer between 10 and 30 and N\* represents additional non-homologous nucleotides which do not base pair with the corresponding genomic sequence in the respective allele and n\* is an integer between 5 and 60, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.
2. A method as claimed in claim 1 and which comprises the detection of further human CFTR gene alleles by, in a separate ARMS reaction, the use of one or more of :  
DF508 mutant primer 5'(N)nACCATT3',  
3849+10kb C>T mutant primer 5'(N)nTACGCA3',  
N1303K mutant primer 5'(N)nTCCATC3',  
1717-1G>A mutant primer 5'(N)nTAATTA3',  
W1282X mutant primer 5'(N)nCAGTCA3', and  
G542X mutant primer 5'(N)nTTCTCT3' wherein N and n are as previously defined, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.
3. A method as claimed in claim 1 and which comprises the detection of further human CFTR gene alleles by, in a separate ARMS reaction, the use of one or more of :  
DF508 non-mutant primer 5'(N)nACCACA3',

- 39 -

W1282X mutant primer 5'(N)nCAGTCA3',

1717-1 mutant primer 5'(N)nTAATTA3',

G542X mutant primer 5'(N)nTTCTCT3',

N1303K mutant primer 5'(N)nTCCATC3',

5 DF508 non-mutant primer 5'(N)nACCACA3',

DF508 mutant primer 5'(N)nACCATT3' and

3849+10kb C&gt;T mutant primer 5'(N)nTACGCA3' wherein N and n are as previously

defined, in the presence of appropriate nucleotide triphosphates and an agent for

polymerisation, such that a diagnostic primer is extended only when the corresponding allelic

10 variant is present in the sample; and detecting the presence or absence of the allelic variant by

reference to the presence or absence of a diagnostic primer extension product.

4. A method as claimed in claim 1 and which comprises the detection of further human CFTR gene alleles by, in a separate ARMS reaction, the use of one or more of:

15 A455E mutant primer 5'(N)nGTTGTA3',

2183AA&gt;G mutant primer 5'(N)nGATAGC3',

3659delC mutant primer 5'(N)nCCTAGA3',

DI507 mutant primer 5'(N)nATAACT3',

1078delT mutant primer 5'(N)nTTCCTG3',

20 R347P mutant primer 5'(N)nTCTACC3',

S1251N mutant primer 5'(N)nGAAGCA3' and

E60X mutant primer 5'(N)nCAGTTA3'

wherein N and n are as previously defined, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only

25 when the corresponding allelic variant is present in the sample; and detecting the presence or

absence of the allelic variant by reference to the presence or absence of a diagnostic primer

extension product.

5. A method as claimed in claim 1 and which comprises the detection of further human

30 CFTR gene alleles by, in a separate ARMS reaction, the use of one or more of:

G85E mutant primer 5'(N)nCTACGA3',

405+1G&gt;A mutant primer 5'(N)nTAGTGA3',

S549R mutant primer 5'(N)nCTGACG3',

W1089X mutant primer 5'(N)nCAAATA3' and

D1152H mutant primer 5'(N)nCACTTG3' wherein N and n are as previously defined,

in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such

5 that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

6. A method as claimed in any one of the previous claims and which further comprises  
10 the use of one or more common amplification primer(s) in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, and subjecting the mixture to PCR amplification such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a corresponding PCR amplification product.

15

7. A method as claimed in any one of claims 2-6 and wherein two or more diagnostic primers are used as a multiplex.

8. A method as claimed in any one of claims 2-6 and wherein all of the diagnostic  
20 primers are used in a single multiplex reaction.

9. A method as claimed in any one of the previous claims and which further comprises the use of one or more control primers.

25 10. A method as claimed in claim 1 and used in conjunction with any known diagnostic CFTR gene procedure.

11. A diagnostic primer selected from any one of

5T variant primer 5'(N)nTGTTAAAGAC3',

30 7T variant primer 5'(N\*)n\*(N)nTTAAAAAAGC3' and

9T variant primer 5'(N\*)n\*(N)nAAAAAAAATC3' wherein N, n, N\* and n\* are as defined in claim 1 and n is an integer between 6 and 26

12. A diagnostic primer selected from any one of  
 5T variant primer 5'TAATCCCCAAATCCCTGTAAAGAC3',  
 7T variant primer 5'(N\*)n\*TAATCCCCAAATCCCTGTAAAAAAGC3' and  
 5 9T variant primer 5'(N\*)n\*TAATCCCCAAATCCCTGTAAAAAAAATC3' wherein N\*  
 and n\* are as defined in claim 1.
13. A diagnostic primer selected from any one of  
 7T variant primer  
 10 5'GTTAATCATTACGCTACTACGCACCTAATCCCCAAATCCCTGTAAAAAAGC3'  
 and 9T variant primer  
 5'GACTGTACGATACTCATTTATATGAAGTCAGCTACTTACCTATAGAACGCTTGC  
 TAGTTTAATCCCCAAATCCCTGTAAAAAAAATC3'
- 15 14. A set of diagnostic primers as set out in Table 1
15. A diagnostic primer selected from any one of  
 DF508 mutant primer 5'(N)nACCATT3',  
 3849+10kb C>T mutant primer 5'(N)nTACGCA3',  
 20 N1303K mutant primer 5'(N)nTCCATC3',  
 1717-1G>A mutant primer 5'(N)nTAATTA3',  
 W1282X mutant primer 5'(N)nCAGTCA3', and  
 G542X mutant primer 5'(N)nTTCTCT3' wherein N and n are as defined in claim 1 and n is  
 an integer between 6 and 26
- 25
16. A diagnostic primer selected from any one of  
 W1282X mutant primer 5'(N)nGCAACAGTCA3',  
 1717-1G>A mutant primer 5'(N)nTTGGTAATTA3',  
 G542X mutant primer 5'(N)nATAGTTCTCT3',  
 30 N1303K mutant primer 5'(N)nGGGATCCATC3',  
 DF508 mutant primer 5'(N)nAAACACCAT3' and  
 3849+10kb C>T mutant primer 5'(N)nGTCTTACGCA3' wherein N and n are as defined in  
 claim 1 and n is an integer between 6 and 26



17. A diagnostic primer selected from any one of  
W1282X mutant primer 5'TCTTGGGATTCAATAACTTTGCAACAGTCA3',  
1717-1G>A mutant primer 5'TCTCGAATTTTCTATTTTGGTAATTA3',  
5 G542X mutant primer 5'AGTTTGCAGAGAAAGACAATATAGTTCTCT3',  
N1303K mutant primer 5'TGATCACTCCACTGTTTCATAGGGATCCATC3',  
DF508 mutant primer 5'GTATCTATATTCATCATAGGAAACACCATT3', and  
3849+10kb C>T mutant primer 5'GAACATTTCCTTTCAGGGTGTCTTACGCA3'.
- 10 18. A set of diagnostic primers as set out in Table 2
19. A diagnostic primer selected from any one of  
DF508 non-mutant primer 5'(N)nACCACA3',  
W1282X mutant primer 5'(N)nCAGTCA3',  
15 1717-1 mutant primer 5'(N)nTAATTA3',  
G542X mutant primer 5'(N)nTTCTCT3',  
N1303K mutant primer 5'(N)nTCCATC3',  
DF508 non-mutant primer 5'(N)nACCACA3',  
DF508 mutant primer 5'(N)nACCATT3' and  
20 3849+10kb C>T mutant primer 5'(N)nTACGCA3' wherein N and n are as defined in claim 1  
and n is an integer between 6 and 26
20. A diagnostic primer selected from any one of  
DF508 non-mutant primer 5'(N)nAAACACCACA3',  
25 R117H mutant primer 5'(N)nGCGATAGACT3',  
621+1G>T mutant primer 5'(N)nGAAGTATTGA3',  
R334W mutant primer 5'(N)nATCATCCTGT3',  
R1162X mutant primer 5'(N)nTCTGTGAGTT3',  
R553X mutant primer 5'(N)nTTCTTGCTGA3' and  
30 G551D mutant primer 5'(N)nGCTCGTTGTT3' wherein N and n are as defined in claim 1  
and n is an integer between 6 and 26

21. A diagnostic primer selected from any one of

R117H mutant primer 5'AGCCTATGCCTAGATAAATCGCGATAGACT3',  
621+1G>T mutant primer 5'TGCCATGGGGCCTGTGCAAGGAAGTATTGA3',  
R334W mutant primer 5'CCTATGCACTAATCAAAGGAATCATCCTGT3',  
5 R1162X mutant primer 5'TATTTTATTTTCAGATGCGATCTGTGAGTT3',  
R553X mutant primer 5'TTATTCACCTTGCTAAAGAAATTCTTGCTGA3',  
G551D mutant primer 5'GCTAAAGAAATTCTTGCTCGTTGTT3'.

22. A set of diagnostic primers as set out in Table 3

10

23. A diagnostic primer selected from any one of

A455E mutant primer 5'(N)nGTTGTA3',  
2183AA>G mutant primer 5'(N)nGATAGC3',  
3659delC mutant primer 5'(N)nCCTAGA3',  
15 DI507 mutant primer 5'(N)nATAACT3',  
1078delT mutant primer 5'(N)nTTCCTG3',  
R347P mutant primer 5'(N)nTCTACC3',  
S1251N mutant primer 5'(N)nGAAGCA3' and  
E60X mutant primer 5'(N)nCAGTTA3' wherein N and n are as defined in claim 1 and n is  
20 an integer between 6 and 26

24. A diagnostic primer selected from any one of

A455E mutant primer 5'(N)nAGTTGTTGTA3',  
1078delT mutant primer 5'(N)nAGGGTTCCTG3',  
25 R347P mutant primer 5'(N)nTTGTTCTACC3',  
DI507 mutant primer 5'(N)nGAAAATAACT3',  
3659delC mutant primer 5'(N)nTAAACCTAGA3',  
2183AA>G mutant primer 5'(N)nAAAAGATAGC3',  
S1251N mutant primer 5'(N)nCAGGGAAGCA3' and

30 E60X mutant primer 5'(N)nAAGCCAGTTA3' wherein N and n are as defined in claim 1 and  
n is an integer between 6 and 26

25. A diagnostic primer selected from any one of  
A455E mutant primer 5'TTCAAGATAGAAAGAGGACAGTTGTTGTA3',  
1078delT mutant primer 5'CCTTCTTCTTCTCAGGGTTCCTG3',  
R347P mutant primer 5'CACCATCTCATTCTGCATTGTTCTACC3',  
5 DI507 mutant primer 5'GCCTGGCACCATTAAAGAAAATAACT3',  
3659delC mutant primer 5'ATGCCAACAGAAGGTAAACCTAGA3',  
2183AA>G mutant primer 5'CAAACCTCTCCAGTCTGTTTAAAAGATAGC3',  
S1251N mutant primer 5'GGAAGAACTGGATCAGGGAAGCA3' and  
E60X mutant primer 5'TTAGGATTTTTCTTTGAAGCCAGTTA3'.

10

26. A set of diagnostic primers as set out in Table 4

27. A diagnostic primer selected from any one of  
G85E mutant primer 5'(N)nCTACGA3',  
15 405+1G>A mutant primer 5'(N)nTAGTGA3',  
S549R mutant primer 5'(N)nCTGACG3',  
W1089X mutant primer 5'(N)nCAAATA3' and  
D1152H mutant primer 5'(N)nCACTTG3' wherein N and n are as defined in claim 1 and n is  
an integer between 6 and 26

20

28. A diagnostic primer selected from any one of  
G85E mutant primer 5'(N)nTGTTCTACGA3',  
405+1G>A mutant primer 5'(N)nTATTTAGTGA3',  
S549R mutant primer 5'(N)nCACTGACG3',  
25 W1089X mutant primer 5'(N)nCTGCCAAATA3',  
D1152H mutant primer 5'(N)nTATCCACTTG3' wherein N and n are as defined in claim 1  
and n is an integer between 6 and 26

29. A diagnostic primer selected from any one of

30 G85E mutant primer  
5'TAGCCATTGATGACGGAGCGATGTTTTTCTGGAGATTTATGTTCTACGA3'  
405+1G>A mutant primer

- 45 -

5'GATTTATGTTCTATGGAATCTTTTATATTTAGTGA3',

S549R mutant primer 5'TGGAGAAGGTGGAATCACACTGACG3',

W1089X mutant primer 5'AAGCTCTGAATTTACATACTGCCAAATA3' and

D1152H mutant primer 5'AAAGATGATAAGACTTACCAAGCTATCCACTTG3'

5

30. A set of diagnostic primers as set out in Table 5

## (12) INTERNATIONAL APPLICATION PUBLISHED UNDER THE PATENT COOPERATION TREATY (PCT)

(19) World Intellectual Property Organization  
International Bureau(43) International Publication Date  
29 March 2001 (29.03.2001)

PCT

(10) International Publication Number  
WO 01/21833 A2

- (51) International Patent Classification<sup>7</sup>: C12Q 1/68 (74) Agent: PHILLIPS, Neil, Godfrey, Alasdair; AstraZeneca, Global Intellectual Property, P.O. Box 272, Mereside, Alderley Park, Macclesfield, Cheshire SK10 4GR (GB).
- (21) International Application Number: PCT/GB00/03597
- (22) International Filing Date:  
19 September 2000 (19.09.2000) (81) Designated States (*national*): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CR, CU, CZ, DE, DK, DM, DZ, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TR, TT, TZ, UA, UG, US, UZ, VN, YU, ZA, ZW.
- (25) Filing Language: English
- (26) Publication Language: English
- (30) Priority Data:  
9922527.8 24 September 1999 (24.09.1999) GB (84) Designated States (*regional*): ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GW, ML, MR, NE, SN, TD, TG).
- (71) Applicant (*for all designated States except MG, US*): ASTRAZENECA AB [SE/SE]; S-151 85 Sodertalje (SE).
- (71) Applicant (*for MG only*): ASTRAZENECA UK LIMITED [GB/GB]; 15 Stanhope Gate, London W1Y 6LN (GB).
- (72) Inventors; and
- (75) Inventors/Applicants (*for US only*): BAYLIFFE, Andrew, Iain [GB/GB]; Alderly Park, Macclesfield, Cheshire SK10 4TG (GB). DOCTER, Eelco [NL/GB]; Alderly Park, Macclesfield, Cheshire SK10 4TG (GB). KELLY, Stephen, James [GB/GB]; Alderly Park, Macclesfield, Cheshire SK10 4TG (GB). ROBERTSON, Nancy, Hastings [GB/GB]; Alderly Park, Macclesfield, Cheshire SK10 4TG (GB).

## Published:

— Without international search report and to be republished upon receipt of that report.

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

(54) Title: ASSAY


(57) Abstract: A diagnostic method for the detection of the 5T, 7T and 9T alleles in intron (8) of the human CFTR gene which method comprises contacting a test sample of nucleic acid from an individual with a multiplex of diagnostic primers comprising (i) 5T variant primer 5'(N)nAAAGAC3', (ii) 7T variant primer 5'(N\*)n\*(N)nAAAAGC3' and (iii) 9T variant primer 5'(N\*)n\*(N)nAAAATC3', wherein N represents additional nucleotides which base pair with the corresponding genomic sequence in the respective allele and n is an integer between 10 and 30 and N\* represents additional non-homologous nucleotides which do not base pair with the corresponding genomic sequence in the respective allele and n\* is an integer between 5 and 60, in the presence of appropriate nucleotide triphosphates and an agent for polymerisation, such that a diagnostic primer is extended only when the corresponding allelic variant is present in the sample; and detecting the presence or absence of the allelic variant by reference to the presence or absence of a diagnostic primer extension product.

WO 01/21833 A2

WO 01/21833

PCT/GB00/03597

- 1 of 5 -

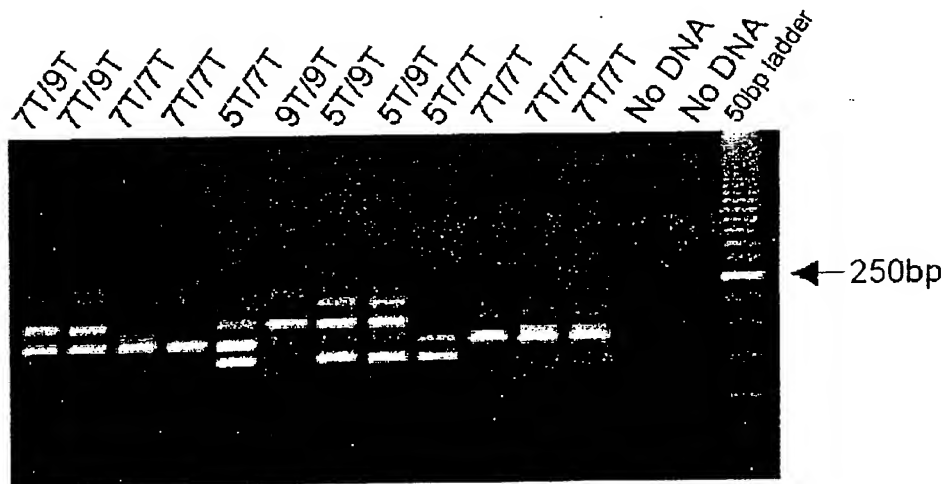


Typical Results from the CF *Poly-T* test

**Figure 1**

5

10

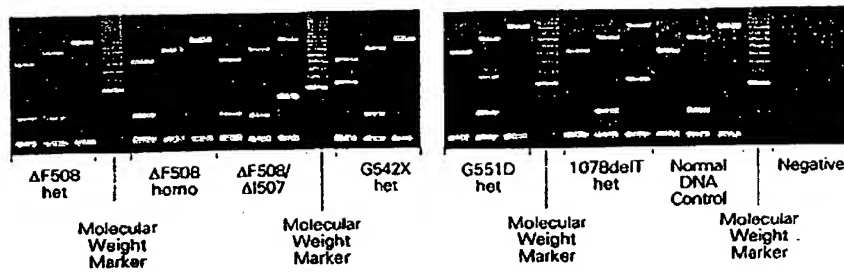

15

**Figure 2**


WO 01/21833

PCT/GB00/03597

- 2 of 5 -

# CF-MEP Agarose Gel Results

D1152H  
W1089X  
G85E  
405+1 G>A  
S549R

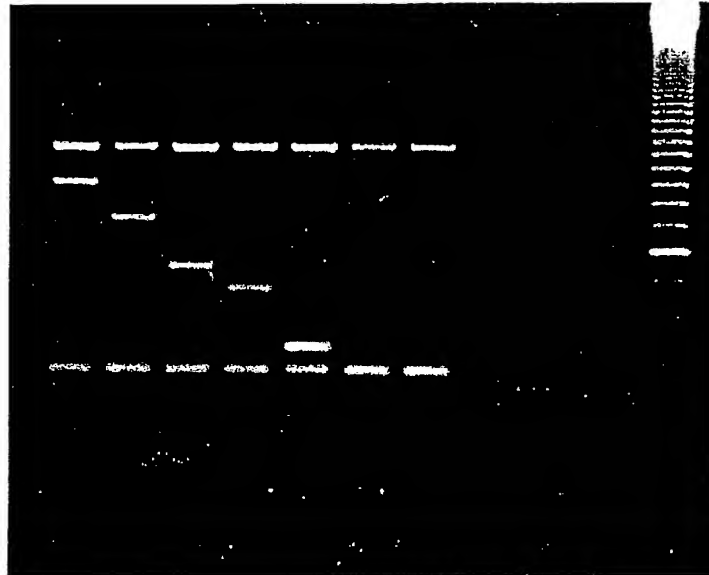


5

Figure 3

WO 01/21833

PCT/GB00/03597

- 3 of 5 -

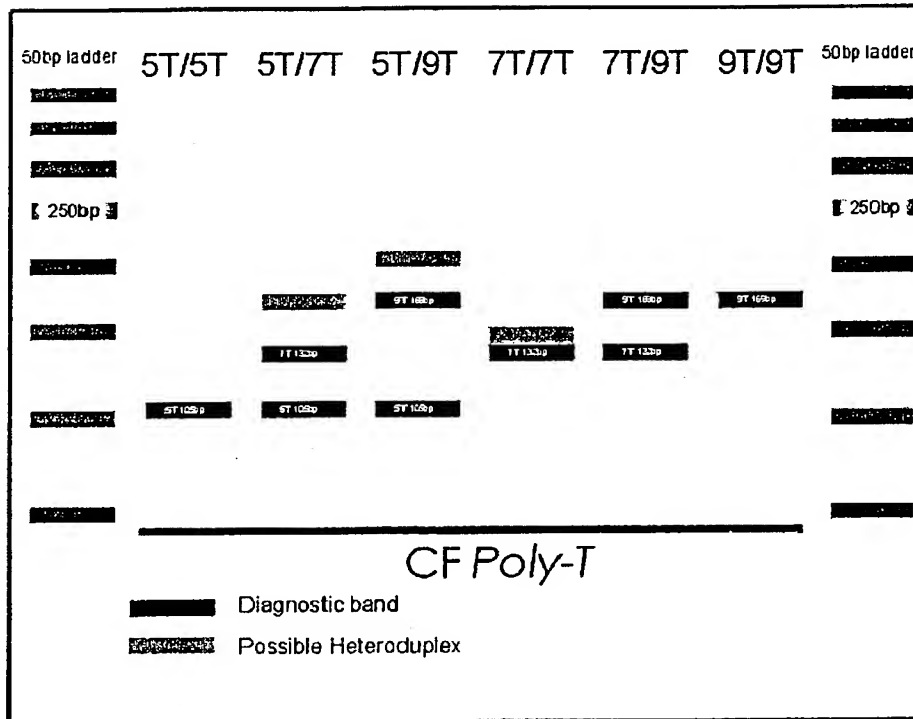


5


**Figure 4**



# ELUCIGENE™ CF20 Test Format

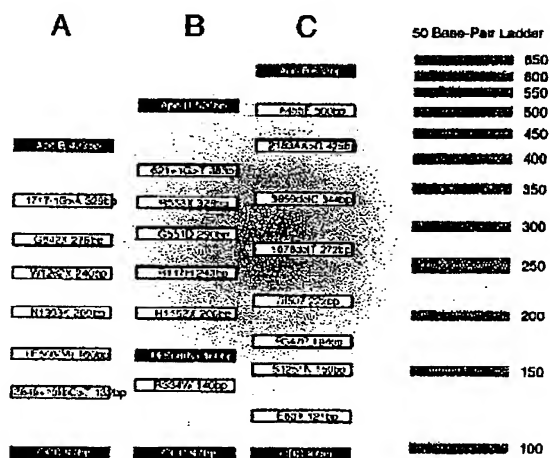

**Figure 5**

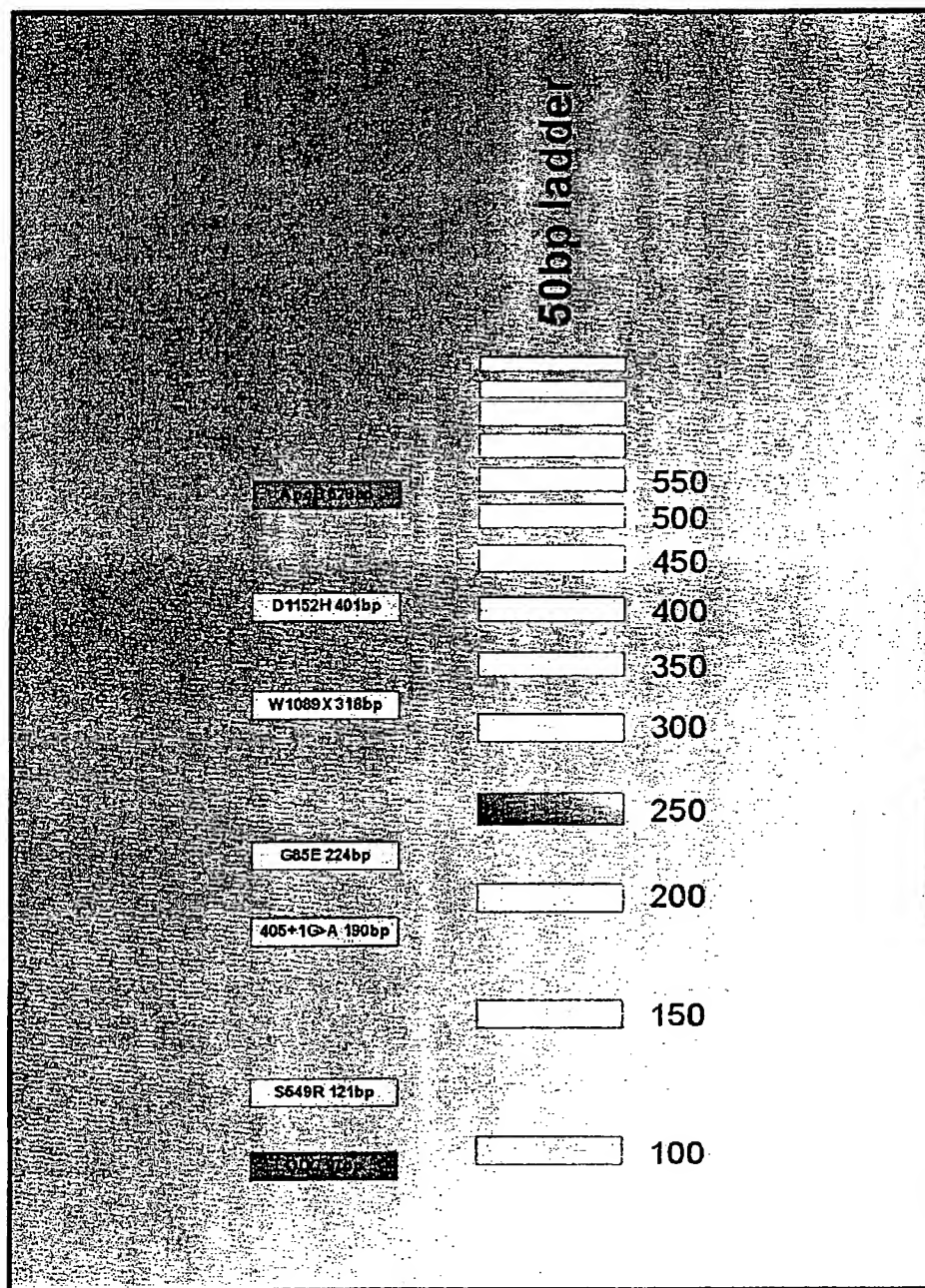

13237.1

13134PCTUS

**DECLARATION AND POWER OF ATTORNEY**

As a below-named inventor, I hereby declare that:

My residence, post office address, and citizenship are as stated below next to my name.

I believe I am the original, first and sole inventor (if only one name is listed below) or an original, first and joint inventor (if plural names are listed below) of the subject matter, which is claimed and for which a patent is sought on the invention entitled:

**'ASSAY'**

the specification of which  
(check one)

☐ is attached hereto

☒ was filed on March 22, 2002 as

Application Serial No. 10/089,001

and was amended on \_\_\_\_\_ (if applicable)

and

☒ described and claimed in PCT International Application No. PCT/GB00/03597 filed on

September 19, 2000.

and was amended on \_\_\_\_\_ under PCT Article 19 (if applicable),

(the undersigned hereby authorizes its attorney to amend this document to insert the filing date and application serial number when they become known.)

I hereby state that I have reviewed and understand the contents of the above-identified specification, including the claims, as amended by any amendment(s) referred to above.

I acknowledge the duty to disclose information of which I am aware which is material to the patentability of this application in accordance with Title 37, Code of Federal Regulations, Section 1.56.

I hereby claim foreign priority benefits under Title 35, United States Code Section 119(a)-(d) or Section 365 of any foreign application(s) for patent, inventor=s certificate, or any PCT international application(s) which designated at least one country other than the United States of America, identified below and have also identified below any foreign application(s) for patent, inventor=s certificate, or any PCT international application designating at least one country other than the United States of America, having a filing date before that of the application(s) on which priority is claimed:

13134PCTUS

## Prior Foreign Application(s)

Number	Country (or PCT)	Filing Date	Priority Claimed	
			Yes	No
9922527.8	GB	24 September 1999	X	

I hereby claim the benefit under Title 35, United States Code, Section 119(e) of any United States provisional application(s) listed below:

Provisional Application No.	Filing Date

I hereby claim the benefit under Title 35, United States Code, Section 120, of any United States application(s), or any PCT international application(s) designating the United States of America, that is/are listed below and, insofar as the subject matter of each of the claims of this application is not disclosed in the prior United States or PCT International application in the manner provided by the first paragraph of Title 35, United States Code, Section 112, I acknowledge the duty to disclose information which is material to patentability, as defined in 37 CFR ' 1.56, which became available between the filing date of the prior application and the national or PCT International filing date of this application.

<u>Application Serial No.</u>	<u>Filing Date</u>	<u>Status</u> (Patented, pending, or abandoned)
<u>PCT/GB00/03597</u>	<u>19 September 2000</u>	<u>pending</u>

And I hereby appoint

Name	Reg. No.	Name	Reg. No.
David A. Kalow	29,397	Milton Springut	27,721
John J. Santalone	32,794	J. David Ellett, Jr.	27,875
William D. Schmidt	39,492	Gary Molnar	30,299
Sylvia Chiou-Tan	47,324	Scott D. Locke	44,877
Tor Smeland	43,131		

each of them of KALOW & SPRINGUT LLP, 488 Madison Avenue, 19<sup>th</sup> floor, New York, New York 10022, my attorneys, each with full power of substitution and revocation, to prosecute this application, to make alterations and amendments therein, to receive the patent, to transact all business in the Patent and Trademark Office connected therewith and to file any International Applications which are based thereon under the provisions of the Patent Cooperation Treaty.

10089001 022003

Please address all communications, and direct all telephone calls, regarding this application to

David A. Kalow  
KALOW & SPRINGUT LLP  
488 Madison Avenue, 19<sup>th</sup> floor  
New York, New York 10022  
(212) 813-1600

I hereby declare that all statements made herein of my own knowledge are true and that all statements made on information and belief are believed to be true; and further that these statements were made with the knowledge that willful false statements and the like so made are punishable by fine or imprisonment, or both, under Section 1001 of Title 18 of the United States Code and that such willful false statements may jeopardize the validity of the application or any patent issued thereon.

Full name of sole or first inventor: Andrew Bayliffe

Citizenship: Great Britain

Residence Address: SK 10 4TG

Alderley Park, Macclesfield, Cheshire, Great Britain *GBN*

Post Office Address: Same As Above

Date of signature 20-2-03

Inventor's signature *A. Bayliffe*

Full name of second inventor: Elco Docter

Citizenship: Netherland

Residence Address: SK 10 4 TG

Alderley Park, Macclesfield, Cheshire, Great Britain *GBN*

Post Office Address: Same As Above

Date of signature 6.01.03

Inventor's signature *[Signature]*

10069001 022003

300  
Full name of third inventor:Stephen James ~~Kelly~~ Kelly

Citizenship:

Great Britain

Residence Address: SK 10 4TGAlderley Park, Macclesfield, Cheshire, Great Britain CBNPost Office Address: Same As AboveDate of signature 6<sup>th</sup> Jan 2003

Inventor's signature

Stephen J. Kelly400  
Full name of fourth inventor:Nancy Hastings Robertson

Citizenship: Great Britain

Residence Address: SK 10 4TGAlderley Park, Macclesfield, Cheshire, Great BritainPost Office Address : Same As AboveDate of signature 6<sup>th</sup> Jan 2003

Inventor's signature

Nancy H. Robertson

WO 01/21833

PCT/GB00/03597

- 1 -

## SEQUENCE LISTING

<110> ASTRAZENECA AB

5 <120> ASSAY

<130> NGAP/PHM70607WO

<140>

10 <141>

<150> 9922527.8

<151> 1999-09-24

15 <160> 85

<170> PatentIn Ver. 2.1

<210> 1

20 <211> 11

<212> DNA

<213> Artificial Sequence

<220>

25 <223> Description of Artificial Sequence: PCR primer

<400> 1

ntgttaaaga c 11

30

<210> 2

<211> 11

<212> DNA

<213> Artificial Sequence

35

<220>

<223> Description of Artificial Sequence: PCR primer

<400> 2

40 nttaaaaaag c 11

<210> 3

<211> 11

45 <212> DNA

<213> Artificial Sequence

<220>

- 2 -

<223> Description of Artificial Sequence: PCR primer

<400> 3  
naaaaaaat c 11

5

<210> 4  
<211> 26  
<212> DNA  
10 <213> Artificial Sequence

<400> 4  
taattcccca aatccctgtt aaagac 26

15

<210> 5  
<211> 29  
<212> DNA  
<213> Artificial Sequence

20

<220>  
<223> Description of Artificial Sequence: PCR primer

<400> 5  
25 ntaattcccc aaatccctgt taaaaaagc 29

<210> 6  
<211> 31  
30 <212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

35

<400> 6  
ntaattcccc aaatccctgt taaaaaaat c 31

40

<210> 7  
<211> 53  
<212> DNA  
<213> Artificial Sequence

45

<220>  
<223> Description of Artificial Sequence: PCR primer

<400> 7



- 3 -

gttaatcatt cagctactac gcacctaatc ccccaaattc ctgttaaaaa agc 53

<210> 8

5 <211> 89

<212> DNA

<213> Artificial Sequence

**<220>**

10 <223> Description of Artificial Sequence: PCR primer

<400> 8

gactgtacga tactcattta tatgaagtca gctacttacc tatagaacgc ttgctagttt 60

aattcccca aatccctgta aaaaaaatc 89

15

<210> 9

<211> 35

<212> DNA

20 <213> Artificial Sequence

**<220>**

<223> Description of Artificial Sequence: PCR primer

25 <400> 9

gtacataaaa caagcatcta ttgaaaatat ctgac 35

<210> 10

30 <211> 11

<212> DNA

&lt;213&gt; Artificial Sequence

**<220>**

35 <223> Description of Artificial Sequence: PCR primer

<400> 10

ngcaacagtc a 11

40

<210> 11

<211> 11

<212> DNA

<213> Artificial Sequence

45

**<220>**

<223> Description of Artificial Sequence: PCR primer

<400> 11  
nttggtatt a

11

```
5  <210> 12
    <211> 11
    <212> DNA
    <213> Artificial Sequence
```

10 <220>  
<223> Description of Artificial Sequence: PCR primer

```
<400> 12
natagttctc t
```

11

15

<210> 13

<211> 11

<212> DNA

20 <213> Artificial Sequence

<223> Description of Artificial Sequence: PCR primer

25 <400> 13  
ngggatccat c

11

```

30      <210> 14
      <211> 11
      <212> DNA
      <213> Artificial Sequence

```

35 <223> Description of Artificial Sequence: PCR primer

```
<400> 14
naaacaccat t
```

11

40

```
<210> 15
<211> 11
<212> DNA
<213> Artificial Sequence
```

45  
 <220>  
 <223> Description of Artificial Sequence: PCR primer

WO 01/21833

PCT/GB00/03597

- 5 -

<400> 15  
ngtcttacgc a 11

5 <210> 16  
<211> 30  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 16  
tcttgggatt caataacttt gcaacagtca 30

15

<210> 17  
<211> 27  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 17  
tctcgaattt tctatttttg gtaatta 27

<210> 18

30 <211> 30  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

35 <400> 18  
agtttgcaga gaaagacaat atagttctct 30

40

<210> 19  
<211> 30  
<212> DNA  
<213> Artificial Sequence

45 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 19  
tgatcactcc actgttcata gggatccatc 30

5 <210> 20  
<211> 30  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 20  
gtatctatat tcatcatagg aaacaccatt 30

15

<210> 21  
<211> 29  
<212> DNA  
20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 21  
gaacatttcc tttcagggtg tcttacgca 29

<210> 22  
30 <211> 11  
<212> DNA  
<213> Artificial Sequence

<220>  
35 <223> Description of Artificial Sequence: PCR primer

<400> 22  
naaacaccac a 11

40

<210> 23  
<211> 11  
<212> DNA  
<213> Artificial Sequence

45

<220>  
<223> Description of Artificial Sequence: PCR primer

WO 01/21833

PCT/GB00/03597

- 7 -

<400> 23  
ngcgatagac t 11

5 <210> 24  
<211> 11  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 24  
ngaagtattg a 11

15

<210> 25  
<211> 11  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 25  
natcatcctg t 11

<210> 26

30 <211> 11  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

35

<400> 26  
ntctgtgagt t 11

40

<210> 27  
<211> 11  
<212> DNA

<213> Artificial Sequence

45

<220>  
<223> Description of Artificial Sequence: PCR primer

WO 01/21833

PCT/GB00/03597

- 8 -

<400> 27  
nttcttgctg a 11

5 <210> 28  
<211> 11  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 28  
ngctcgttgt t 11

15

<210> 29  
<211> 30  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 29  
agcctatgcc tagataaatc gcgatatgact 30

<210> 30

30 <211> 30  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

35 <400> 30  
tgccatgggg cctgtgcaag gaagtattga 30

40

<210> 31  
<211> 30  
<212> DNA

45 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

<400> 31  
cctatgcact aatcaaagga atcatcctgt 30

5 <210> 32  
<211> 30  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 32  
tatttttatt tcagatgcga tctgtgagtt 30

15

<210> 33  
<211> 31  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 33  
ttattcacct tgctaaagaa attcttgctg a 31

<210> 34

30 <211> 25  
<212> DNA  
<213> Artificial Sequence

<220>

35 <223> Description of Artificial Sequence: PCR primer

<400> 34  
gctaaagaaa ttcttgctcg ttgtt 25

40

<210> 35  
<211> 11  
<212> DNA  
<213> Artificial Sequence

45 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 35  
nagttgttgt a

11

5 <210> 36  
<211> 11  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 36  
nagggttcct g

11

15  
  
<210> 37  
<211> 11  
<212> DNA  
20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 37  
nttgttctac c

11

<210> 38  
30 <211> 11  
<212> DNA  
<213> Artificial Sequence

<220>  
35 <223> Description of Artificial Sequence: PCR primer

<400> 38  
ngaaaataac t

11

40  
  
<210> 39  
<211> 11  
<212> DNA  
<213> Artificial Sequence

45  
  
<220>  
<223> Description of Artificial Sequence: PCR primer



<400> 39  
ntaaacctag a

11

```
5  <210> 40
    <211> 11
    <212> DNA
    <213> Artificial Sequence
```

10 <220>  
<223> Description of Artificial Sequence: PCR primer

```
<400> 40
naaaagatag c
```

11

15

<210> 41

<211> 11

<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25      <400> 41  
ncaggggaagc a 11

30      <210> 42  
         <211> 11  
         <212> DNA  
         <213> Artificial Sequence

35 <223> Description of Artificial Sequence: PCR primer

<400> 42  
naagccagtt a

11

40

```
<210> 43
<211> 29
<212> DNA
<213> Artificial Sequence
```

45

<220>

<223> Description of Artificial Sequence: PCR primer

<400> 43  
ttcaagatag aaagaggaca gttgttgta 29

5 <210> 44  
<211> 23  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 44  
cctttcttctt ctcagggttc ctg 23

15

<210> 45  
<211> 27  
<212> DNA  
20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 45  
caccatctca ttctgcattg ttctacc 27

<210> 46

30 <211> 26  
<212> DNA  
<213> Artificial Sequence

<220>  
35 <223> Description of Artificial Sequence: PCR primer

<400> 46  
gcctggcacc attaaagaaa ataact 26

40

<210> 47  
<211> 24  
<212> DNA  
<213> Artificial Sequence

45

<220>  
<223> Description of Artificial Sequence: PCR primer

<400> 47  
atgccaacag aaggtaaacc taga 24

5 <210> 48  
<211> 29  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 48  
caaactctcc agtctgttta aaagatagc 29

15

<210> 49  
<211> 23  
<212> DNA  
20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 49  
ggaagaactg gatcagggaa gca 23

<210> 50  
30 <211> 26  
<212> DNA  
<213> Artificial Sequence

<220>  
35 <223> Description of Artificial Sequence: PCR primer

<400> 50  
ttaggatttt tctttgaagc cagtta 26

40

<210> 51  
<211> 11  
<212> DNA  
<213> Artificial Sequence

45 <220>  
<223> Description of Artificial Sequence: PCR primer

WO 01/21833

PCT/GB00/03597

- 14 -

<400> 51  
ntgttctaag a 11

5 <210> 52  
<211> 11  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 52  
ntatttagtg a 11

15

<210> 53  
<211> 11  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 53  
ncacactgac g 11

<210> 54

30 <211> 11  
<212> DNA  
<213> Artificial Sequence

<220>

35 <223> Description of Artificial Sequence: PCR primer

<400> 54  
nctgccaaat a 11

40

<210> 55  
<211> 11  
<212> DNA  
<213> Artificial Sequence

45 <220>  
<223> Description of Artificial Sequence: PCR primer

WO 01/21833

PCT/GB00/03597

- 15 -

<400> 55  
ntatccactt g

11

5 <210> 56  
<211> 50  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 56  
tagccattga tgacggagcg atgttttttc tggagattta tgttctacga

50

15  
  
20 <210> 57  
<211> 36  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 57  
gatttatggt ctatggaatc tttttatatt tagtga

36

30 <210> 58  
<211> 25  
<212> DNA  
<213> Artificial Sequence

35 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 58  
tggagaaggt ggaatcacac tgacg

25

40  
  
45 <210> 59  
<211> 28  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

<400> 59  
aagctctgaa ttacatact gccaaata 28

5 <210> 60  
<211> 33  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 60  
aaagatgata agacttacca agctatccac ttg 33

15

<210> 61  
<211> 29  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 61  
taatctctac caaatctgga tactatacc 29

<210> 62

30 <211> 29  
<212> DNA  
<213> Artificial Sequence

<220>

35 <223> Description of Artificial Sequence: PCR primer

<400> 62  
taatctctac caaatctgga tactatacc 29

40

<210> 63  
<211> 25  
<212> DNA

<213> Artificial Sequence

45

<220>  
<223> Description of Artificial Sequence: PCR primer

<400> 63  
gaattcccaa acttttagag acatc 25

5 <210> 64  
<211> 30  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 64  
cttgatggta agtacatggg tttttcttat 30

15

<210> 65  
<211> 30  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 65  
ccagacttca cttctaatta tgattatggg 30

<210> 66

30 <211> 31  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

35 <400> 66  
ttgtggatca aatttcagtt gacttgcat c 31

40

<210> 67  
<211> 30  
<212> DNA

<213> Artificial Sequence

45 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 67  
gacttcactt ctaatgatga ttatgggaga 30

5 <210> 68  
<211> 37  
<212> DNA  
<213> Artificial Sequence

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 68  
gtttcacata gtgtatgacc ctctatatac actcatt 37

15

<210> 69  
<211> 31  
<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 69  
atctaaaatt ggagcaatgt tgtttttgac c 31

<210> 70

30 <211> 31  
<212> DNA  
<213> Artificial Sequence

<220>

35 <223> Description of Artificial Sequence: PCR primer

<400> 70  
atctaaaatt ggagcaatgt tgtttttgac c 31

40

<210> 71  
<211> 37  
<212> DNA  
<213> Artificial Sequence

45

<220>  
<223> Description of Artificial Sequence: PCR primer



<400> 71  
gtttcacata gtgtatgacc ctctatatac actcatt

37

```
5  <210> 72
    <211> 29
    <212> DNA
    <213> Artificial Sequence
```

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 72  
ttttgctgtg agatctttga cagtcattt

29

15

<210> 73

<211> 30

<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 73  
tttgtttatt gctccaagag agtcatacca

30

```

30      <210> 74
      <211> 49
      <212> DNA
      <213> Artificial Sequence

```

35 <223> Description of Artificial Sequence: PCR primer

<400> 74  
gactgactga ctgactgaaa tggagacttt ttgtttatgt ggttactaa

49

40

<210> 75  
<211> 35  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

- <400> 75  
gtatgataga gattatatgc aataaaacat taaca 35
- 5 <210> 76  
<211> 34  
<212> DNA  
<213> Artificial Sequence
- 10 <220>  
<223> Description of Artificial Sequence: PCR primer
- <400> 76  
tgtgtctaatt attgattcta ctgtacaata ataa 34
- 15
- <210> 77  
<211> 35  
<212> DNA
- 20 <213> Artificial Sequence
- <220>  
<223> Description of Artificial Sequence: PCR primer
- 25 <400> 77  
atTTTTccaa acttcattag aactgatcta ttgac 35
- <210> 78
- 30 <211> 27  
<212> DNA  
<213> Artificial Sequence
- <220>
- 35 <223> Description of Artificial Sequence: PCR primer
- <400> 78  
cacagtagct tacccataga ggaaaca 27
- 40
- <210> 79  
<211> 35  
<212> DNA
- <213> Artificial Sequence
- 45
- <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 79  
atttttccaa acttcattag aactgatcta ttgac

35

```
5  <210> 80
    <211> 23
    <212> DNA
    <213> Artificial Sequence
```

10 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 80  
gctcacctgt ggtatcactc caa

23

15

<210> 81

<211> 34

<212> DNA

20 <213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

25 <400> 81  
aatcaaaacta tgttaaggga aataggacaa ctaa

34

30      <210> 82  
         <211> 39  
         <212> DNA  
         <213> Artificial Sequence

35 <220>  
<223> Description of Artificial Sequence: PCR primer

<400> 82  
cgattcgatt cagtttctg tggttctta gtgttgga

39


40

<210> 83  
<211> 33  
<212> DNA  
<213> Artificial Sequence

<220>  
<223> Description of Artificial Sequence: PCR primer

- <400> 83  
gtaatttttt tacatgaatg acatttacag caa 33
- 5 <210> 84  
<211> 33  
<212> DNA  
<213> Artificial Sequence
- 10 <220>  
<223> Description of Artificial Sequence: PCR primer
- <400> 84  
ggaaattatt tgtttaacaa taaaacaatg gaa 33
- 15
- <210> 85  
<211> 26  
<212> DNA
- 20 <213> Artificial Sequence
- <220>  
<223> Description of Artificial Sequence: PCR primer
- 25 <400> 85  
ccaacaacac ctccaatacc agtaac 26

United States Patent & Trademark Office  
Office of Initial Patent Examination -- Scanning Division

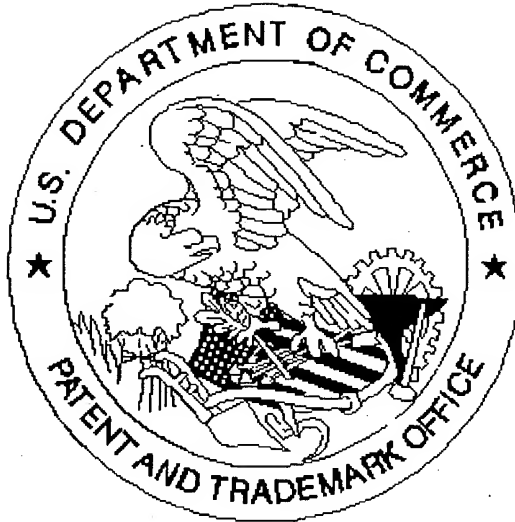

Application deficiencies found during scanning:

☐ Page(s) \_\_\_\_\_ of \_\_\_\_\_ were not present  
for scanning. (Document title)

☐ Page(s) \_\_\_\_\_ of \_\_\_\_\_ were not  
present  
for scanning. (Document title)

*X Scanned copy is best available. Some drawings are dark*